

**Reagent Cover Position Adjustment**

**Reagent Cover Position Adjustment**..... 1

1. Information..... 2

2. Adjustment tool..... 2

3. Preparation..... 2

    3.1 Preparing the adjustment..... 2

    3.2 Start SSW ..... 4

4. Adjustment of Reagent Pipetting Cover A position ..... 4

5. Adjustment of Reagent Pipetting Cover B position ..... 9

6. Accomplishing the adjustment..... 9

1. Information

Make sure the following adjustments.

- Adjusting Reagent Probe R1-A Horizontal (Refer to “S03-03 Reagent Probe Position Adjustment”)
  - Adjusting Reagent Probe R1-B Horizontal (Refer to “S03-03 Reagent Probe Position Adjustment”)
  - Adjusting Reagent Probe R2/3-A Horizontal (Refer to “S03-03 Reagent Probe Position Adjustment”)
  - Adjusting Reagent Probe R2/3-B Horizontal (Refer to “S03-03 Reagent Probe Position Adjustment”)
  - Adjustment of Reagent Disk A (Refer to “S03-08 Reagent Manager Position and Gripper Assy Adjustment”)
  - Adjustment of Reagent Disk B (Refer to “S03-08 Reagent Manager Position and Gripper Assy Adjustment”)
- If they haven't been adjusted yet, perform the adjustments.

2. Adjustment tool


The following adjustment tools are used for this adjustment. Please refer to "S03-01 List of adjustment tools" for a list of jigs.

| P/N      | Name                       | Appearance                                                                            |
|----------|----------------------------|---------------------------------------------------------------------------------------|
| 848-3884 | JACKET ADJUSTING TOOL c703 |  X2 |

3. Preparation

3.1 Preparing the adjustment.

- (1) Open COVER TOP F ASSY (a).

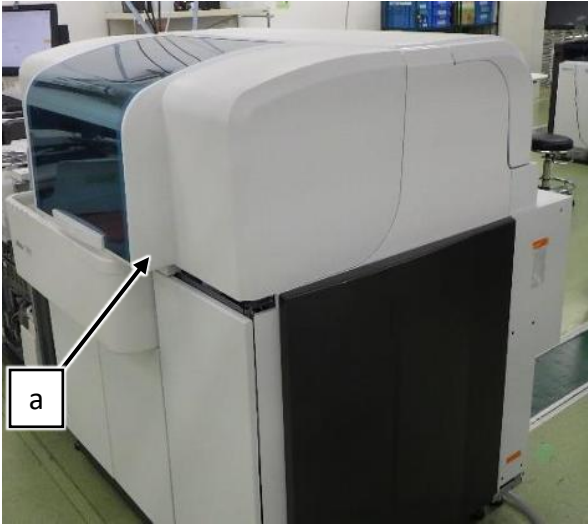


- (2) Loosen the screw and remove COVER MSW ASSY (a).




- (3) Turn on Maintenance switch.



- (4) Remove the 4 screw caps of COVER JACKET A 703 ASSY (a).

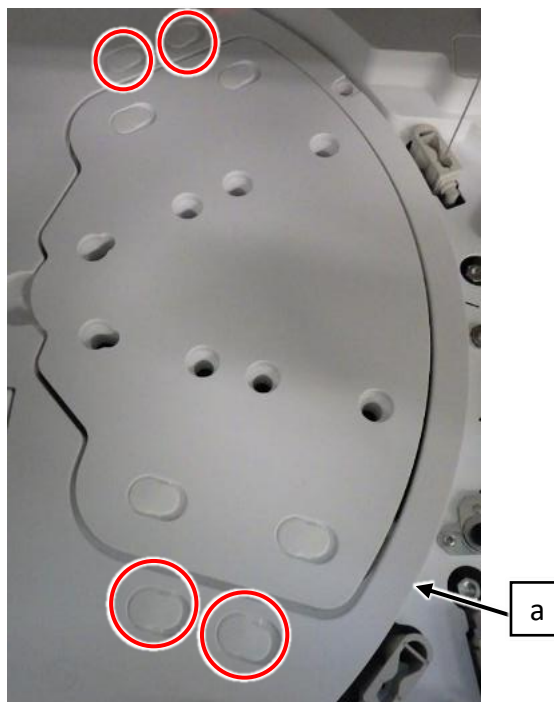


- (5) Remove the 4 screw caps of COVER JACKET B 703 ASSY (a).

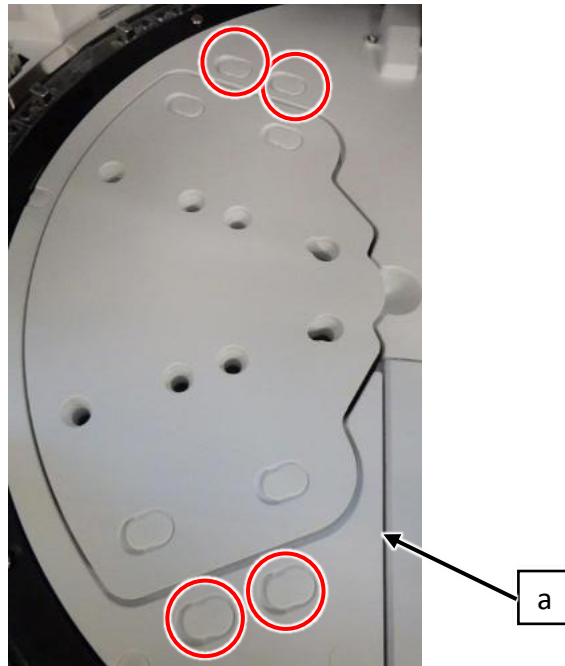

### 3.2 Start SSW


- (1) Refer to “17 pro v2 SSW Adjustment Procedure”.
- (2) Select [Service] > [Adjustment] > [Perform...] in SSW.

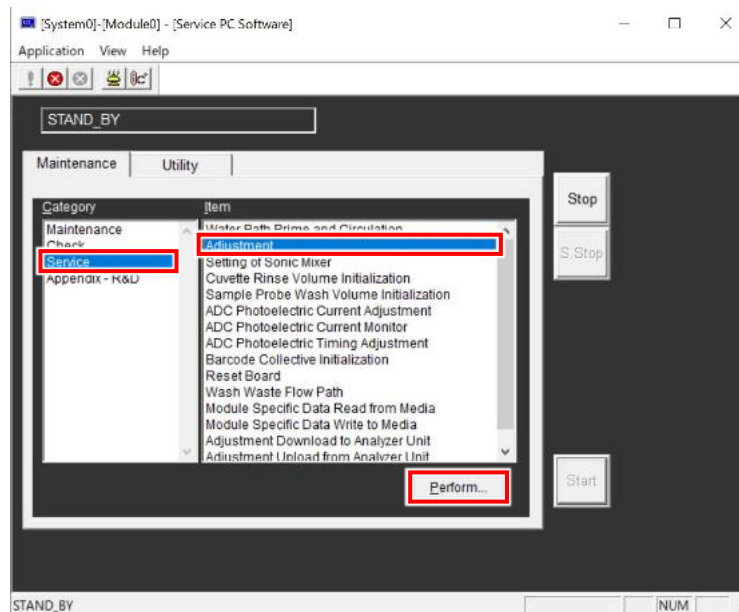


### 4. Adjustment of Reagent Pipetting Cover A position

The adjustment procedure for Reagent Pipetting Cover B is basically same as Reagent Pipetting Cover A.

- (1) Select [Reagent Pipetting cover A] from Adjustment item and select [Adjust...].

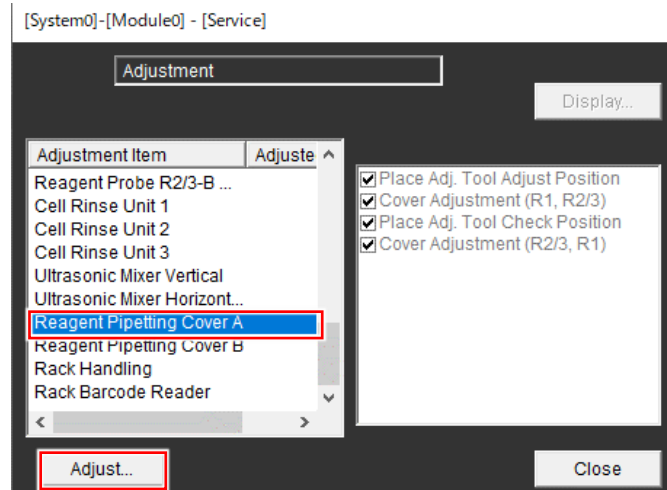


- (2) Set the 2 tools (a) on R1A/position 1 and R23A/position 6 of REAGENT PIPETTING COVER A (b).

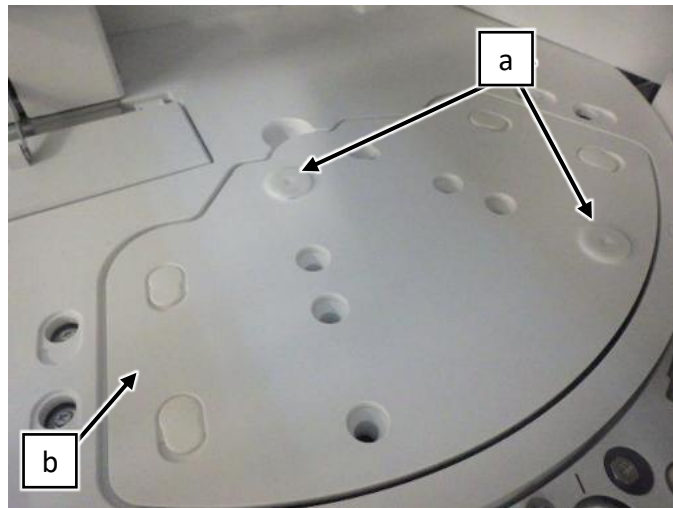


(a): JACKET ADJUSTING TOOL c703 (P/N: 848-3884)


- (3) Select [Set Up].

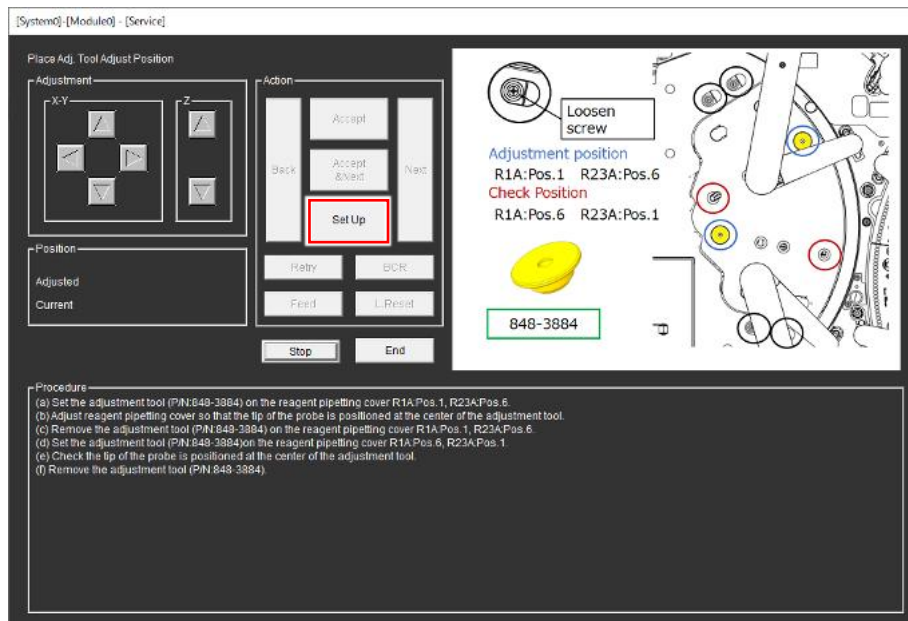


- (4) Lower 2 Reagent probe arms and make sure that each tip of Reagent probe is at the center of the tool.  
If each tip is at the center, proceed to the next step.  
If not, proceed to the step (8).

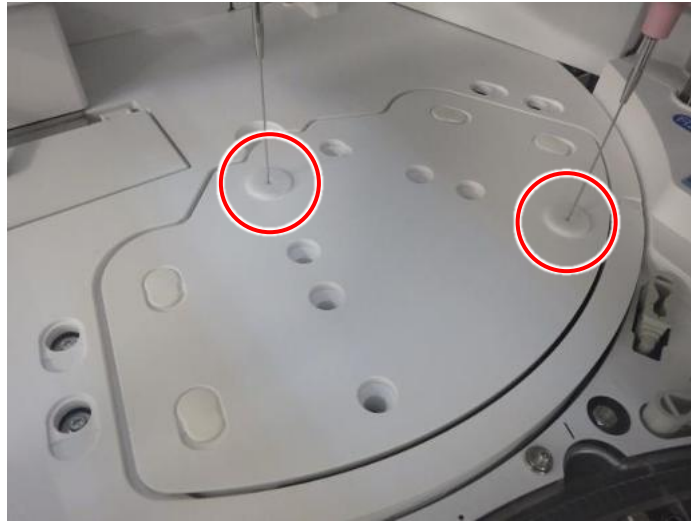

- (5) Loosen the 4 screws.

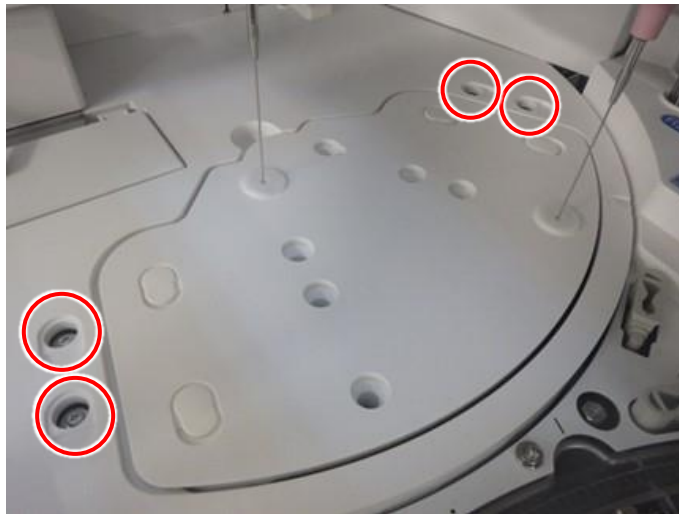

- (6) Adjust REAGENT PIPETTING COVER A (a) position so that each tip of Reagent probe is at the center of the tool.
- (7) Tighten the 4 adjustment screws.

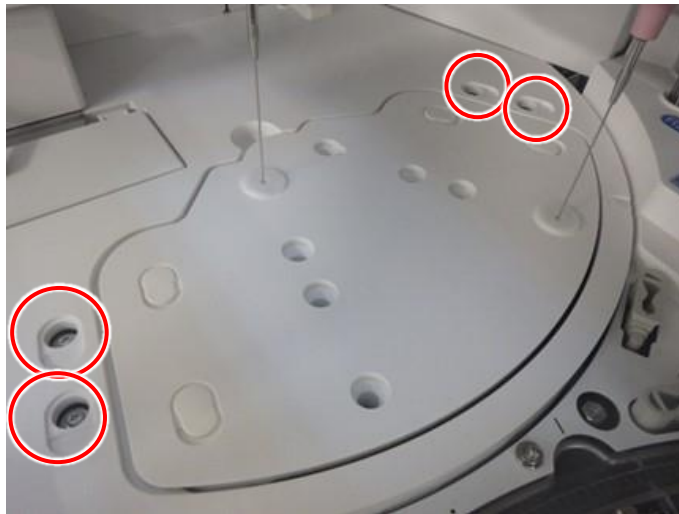

- (8) Lift 2 Reagent probe arms.
- (9) Remove the 2 tools.
- (10) Select [Next].

[System] - [Module] - [Service]

Cover Adjustment (R1, R2/3)

Adjustment

X-Y

Z

Position

Adjusted

Current

Action

Accept

Accept & Feed

Set Up

Back

Feed

Stop

End

Next

Next

Next

Loosen screw

Adjustment position

R1A: Pos.1 R23A: Pos.6

Check Position

R1A: Pos.6 R23A: Pos.1

848-3884

Procedure

(a) Set the adjustment tool (PN: 848-3884) on the reagent pipetting cover R1A: Pos.1, R23A: Pos.6.

(b) Adjust reagent pipetting cover so that the tip of the probe is positioned at the center of the adjustment tool.


(c) Remove the adjustment tool (PN: 848-3884) on the reagent pipetting cover R1A: Pos.1, R23A: Pos.6.

(d) Set the adjustment tool (PN: 848-3884) on the reagent pipetting cover R1A: Pos.6, R23A: Pos.1.

(e) Check the tip of the probe is positioned at the center of the adjustment tool.


(f) Remove the adjustment tool (PN: 848-3884).

(11) Set the 2 tools (a) on R23A/position 1 and R1A/position 6 of REAGENT PIPETTING COVER A (b).

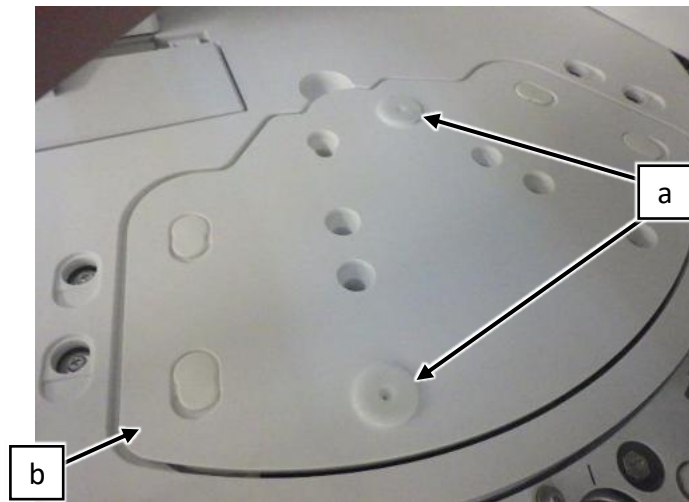

(a): JACKET ADJUSTING TOOL c703 (P/N: 848-3884)

(12) Select [Set Up].

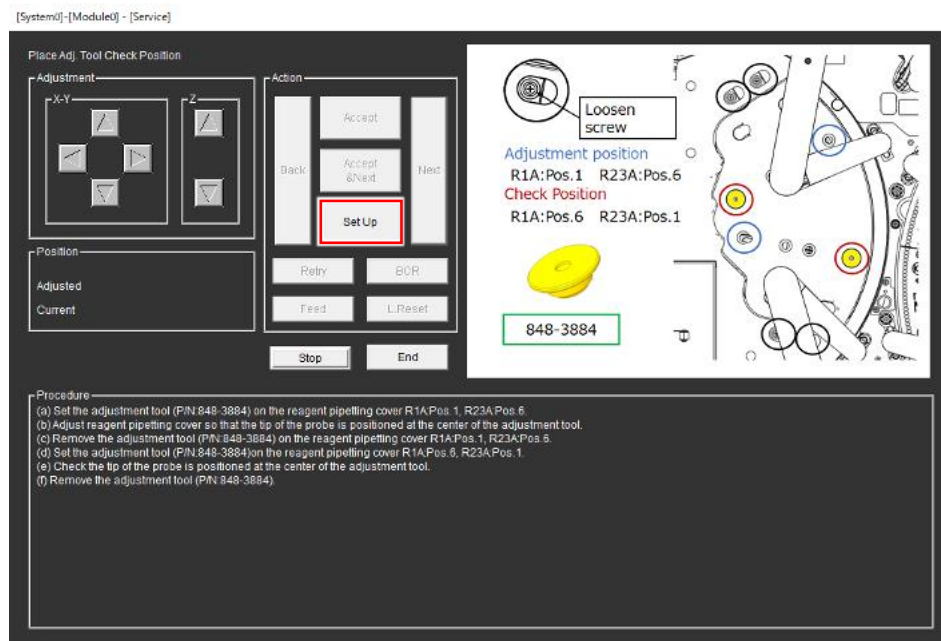




(13) Lower 2 Reagent probe arms and make sure that each tip of Reagent probe is at the center of the tool.

If each tip is at the center, remove the 2 tools and select [End].

If not, remove the 2 tools and suspend this adjustment, and then perform the adjustment of “S03-03 Reagent Probe Position Adjustment”

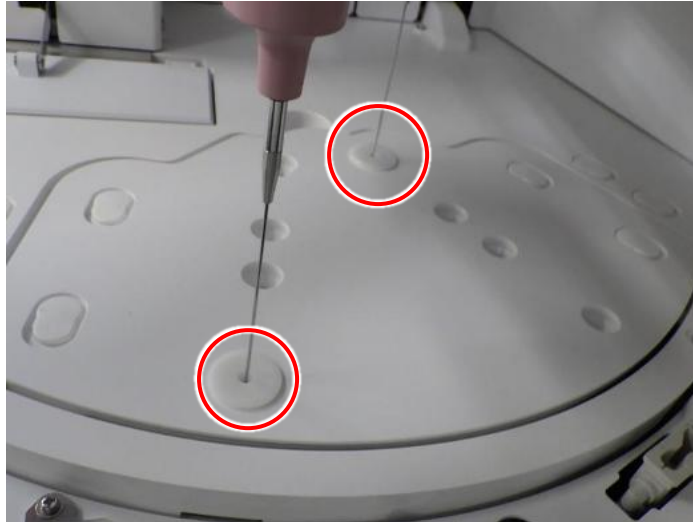

## 5. Adjustment of Reagent Pipetting Cover B position

The adjustment procedure for Reagent Pipetting Cover B is basically same as the procedure for Reagent Pipetting Cover A. Perform Reagent Pipetting Cover B adjustment in the same steps (refer to the chapter 3).

## 6. Accomplishing the adjustment

Refer to “S03-20 Preparing and cleaning up adjustments”.

- (1) Attach the 8 screw caps.
- (2) Turn off Maintenance switch.
- (3) Attach COVER MSW ASSY.
- (4) Close COVER TOP F ASSY.

End of Document