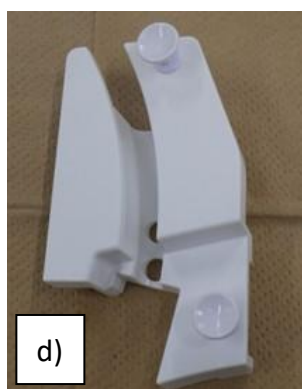
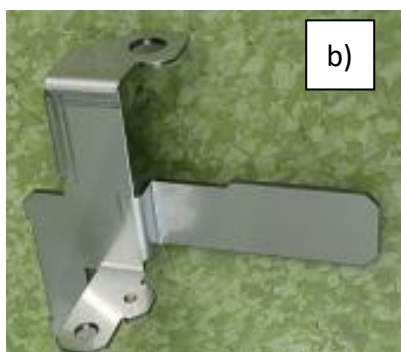
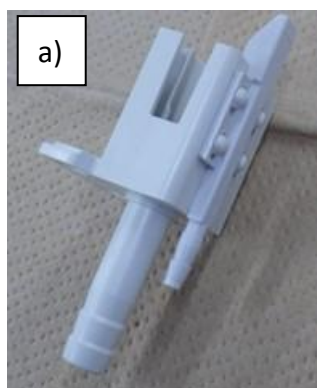


Modification procedure to S1RS MOD (G) KIT (P/N 792-7753)

Applicable Serial Number : up to 2541-10 (Global Ver.)

1. Contents of P/N 792-7753 S1RS MOD (G) KIT

No.	P/N	Parts Name	Quantity	Note
a	792-1324	CLEANING BATH S1 A	1	-
b	792-2777	WS FIXING PLATE 3	1	-
c	792-2881	COVER TB5A	1	-
d	792-0860	COVER TB 7 ASSY	1	-
e	792-7754	TUBE NV2S1RS ASSY	1	-
f	792-7756	JOINT FIX PLATE ASSY	1	-
g	M641135	PAN HEAD SCREW M4X8 ZN3 WITH WASHER	1	-
h	J839050	WIRE BINDER TY-23M	4	-
i	J839051	WIRE BINDER TY-24M	1	-
j	792-1351	SPACER 10	1	-
k	M644039	SCREW M3X10 SUSXM7	2	-



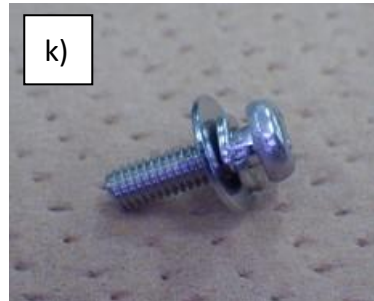


Fig. 1 Contents of P/N 792-7753 S1RS MOD (G) KIT

2. Replacement procedure

- (1) Power off the instruments.

Be sure to wear gloves and eye protection with side shields during work.

- (2) Remove the two nuts to remove the COVER RS CELL.

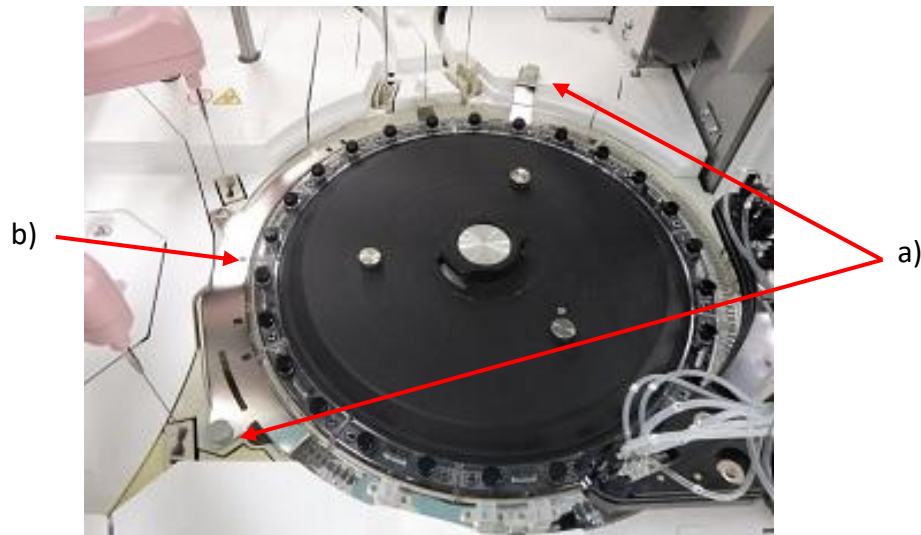


Fig.2 COVER RS CELL

- a) NUT
- b) COVER RS CELL

- (3) Remove the COVER TB 4 ASSY, COVER TB 5 ASSY and COVER TB 7 (4 screws).

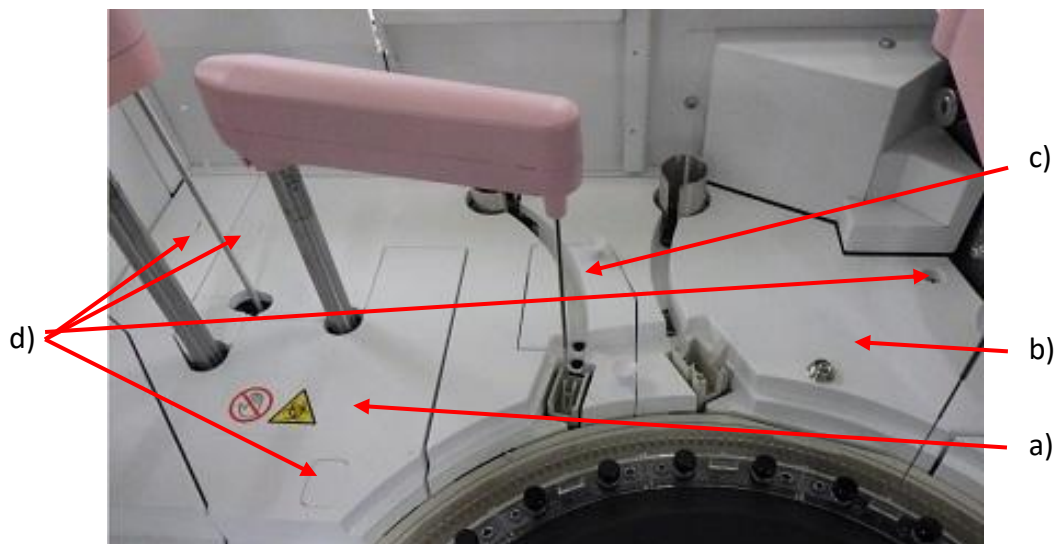


Fig.3 TB COVER

- a) COVER TB 4 ASSY
- b) COVER TB 5 ASSY (remove and keep screws for the new TB COVER 5)
- c) COVER TB 7
- d) Screw

- (4) Remove the COVER HAND RA ASSY and COVER REAR ASSY (8 screws, 1 connector).



Fig.4 COVER HAND RA ASSY

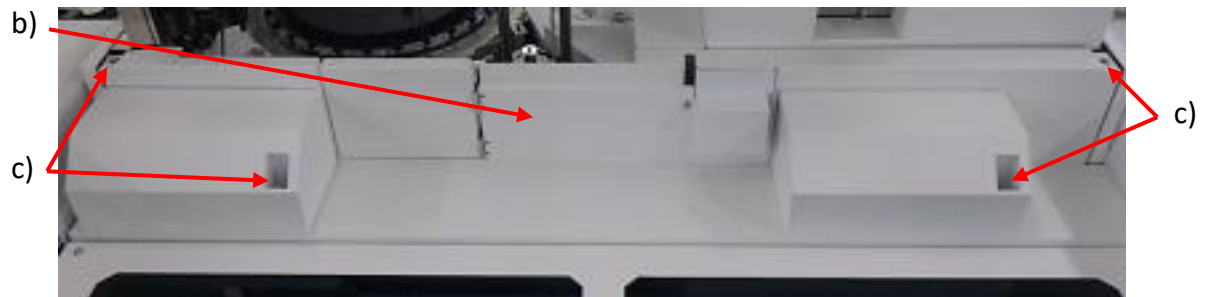


Fig.5 COVER REAR ASSY

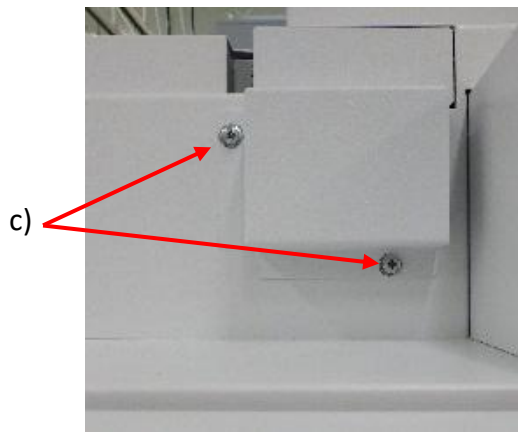


Fig. 6 INTERROCK COVER

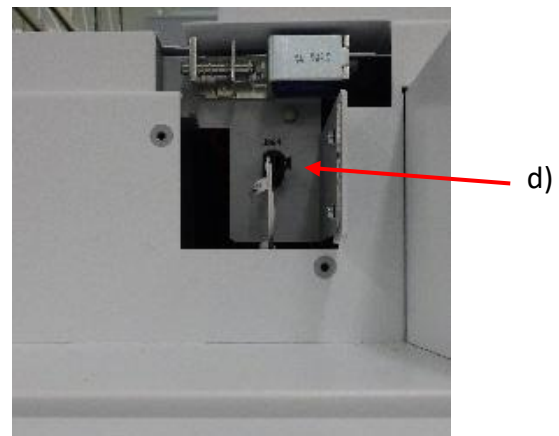


Fig.7 Connector

- a) COVER HAND RA ASSY
- b) COVER REAR ASSY
- c) Screw
- d) Connector

- (5) Remove the Sample Probe 1 for safety
- (6) Remove the COVER SWS 2.(2 screws)

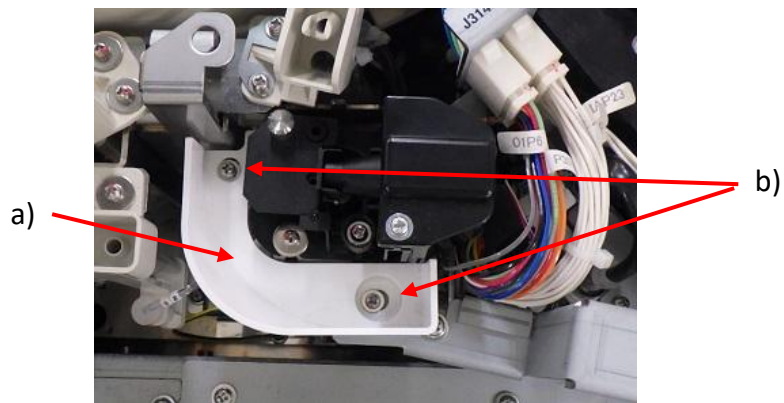


Fig.8 COVER SWS 2

- a) COVER SWS 2
- b) Screw

- (7) Remove the TB COVER FIXING PLATE 8. (2 screws)

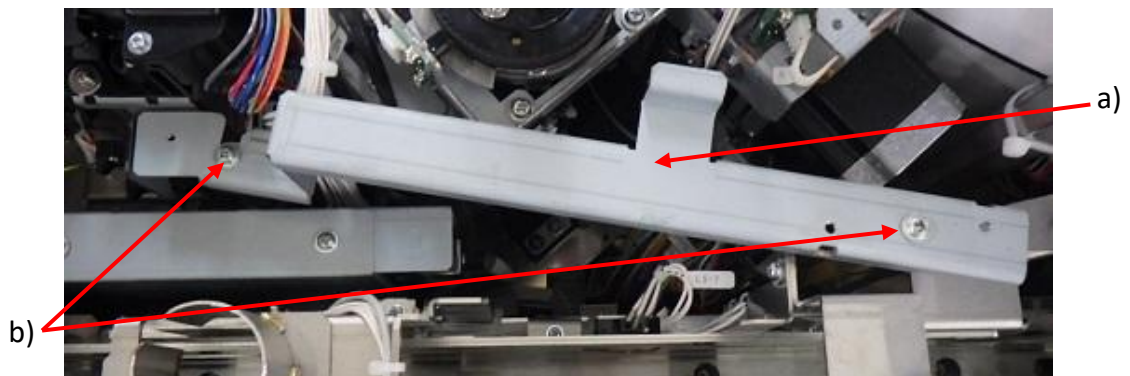


Fig.9 TB COVER FIXING PLATE 8

- a) TB COVER FIXING PLATE 8
- b) Screw

(8) Disconnect the 3 connectors.

Remove the nipple. When a little water flows out of the nipple in removing, wipe it off.

Loosen the 4 screws.

Pull out the PIPETTER S1 ASSY in the upper direction.

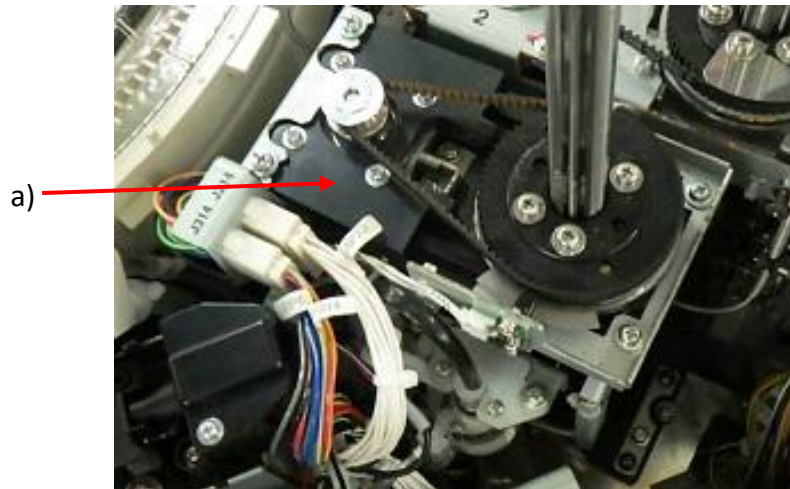


Fig.10 PIPETTER S1 ASSY

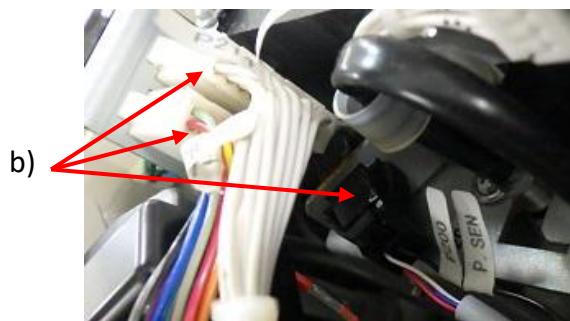


Fig.11 Connector



Fig. 12 Nipple

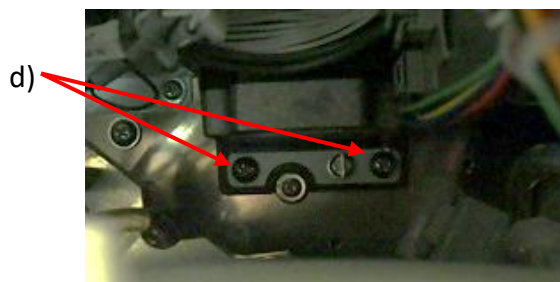


Fig. 13 Screws

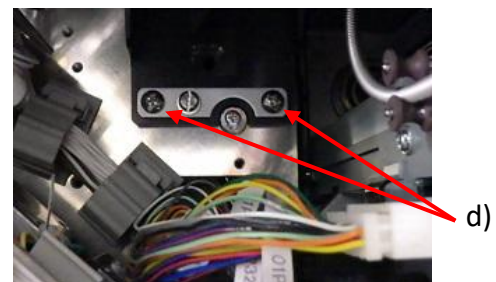


Fig. 14 Screws

- a) PIPETTER S1 ASSY
- b) Connector(J314/J214/J200)
- c) Nipple
- d) Screw

- (9) Remove the two fixing screws of the CLEANING BATH S1.
Remove the SPACER. Discard the removed screws.

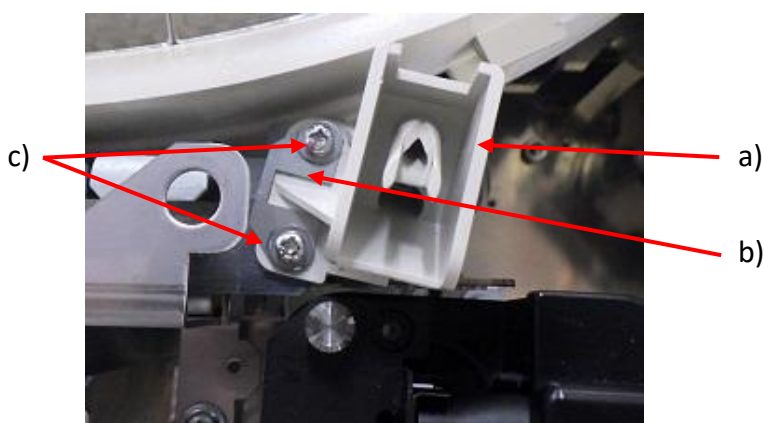


Fig. 15 CLAENING BATH S1

- a) CLAENING BATH S1
- b) SPACER
- c) Screws (Discard the removed screws.)

- (10) Cut the wire binder to remove CLAENING BATH S1 from supply and drain tubes.
Be careful not to damage the tubes by contact with the edges of the metal plate.
Be careful not to damage the CLAENING BATH S1 when removing tube.
Cut a slit vertically in the wire binder head of the tie wrap with cutting tool and pinch the head to remove the wire binder.

Be careful not to damage the drain tube and the surrounding parts with cutting tool.

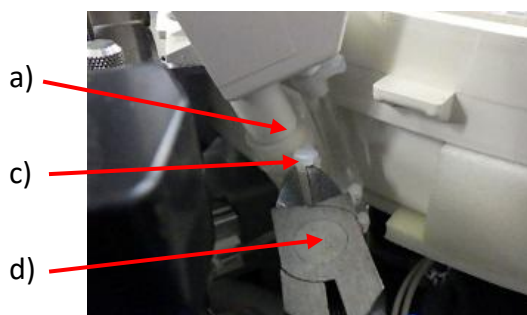


Fig.16 Drain tube

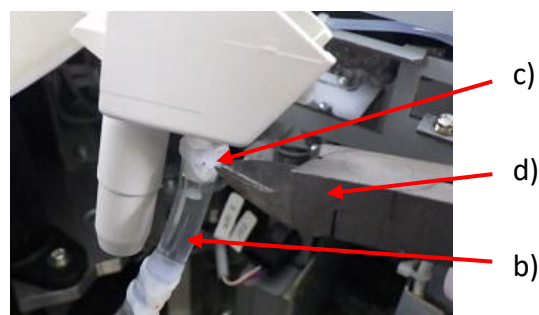


Fig.17 Supply tube

- a) Drain tube
- b) Supply tube
- c) Wire binder head
- d) Cutting tool

(11) Remove the WS FIXING PLATE 2. (1 screw)

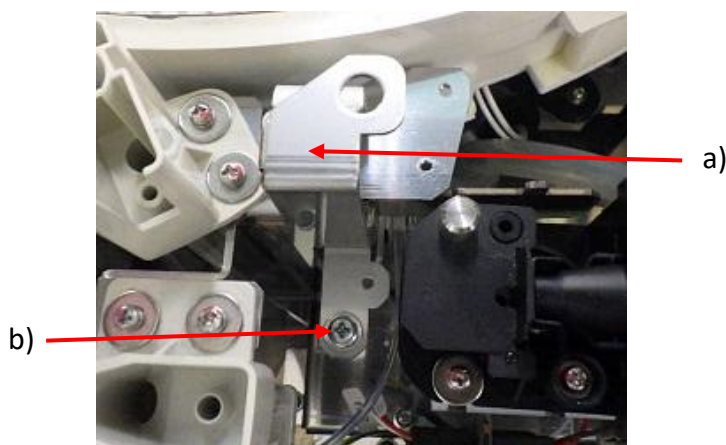


Fig.18 WS FIXING PLATE 2

a) WS FIXING PLATE 2

b) Screw

(12) Insert the supply and drain tubes into the CLEANING BATH S1 A(P/N 792-1324).

The drain tube inserts 15mm, the supply tube inserts to the root.

Secure the connection part with wire binder.

When inserting, align the orientation of CLEANING BATH S1 A (replacement part) with the position of the fixing metal plate to prevent the tubes from being twisted.

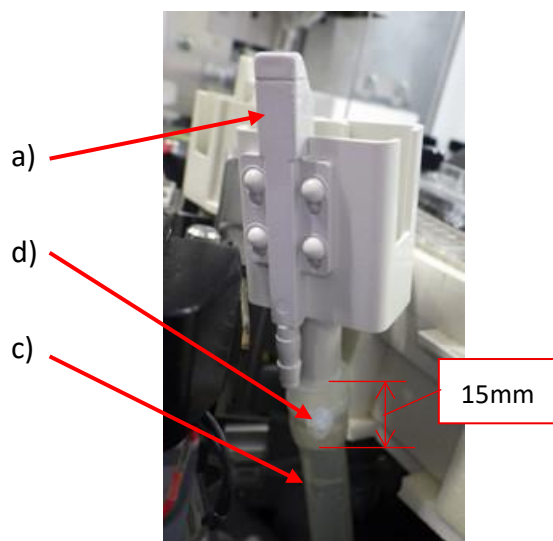


Fig.19 Drain tube

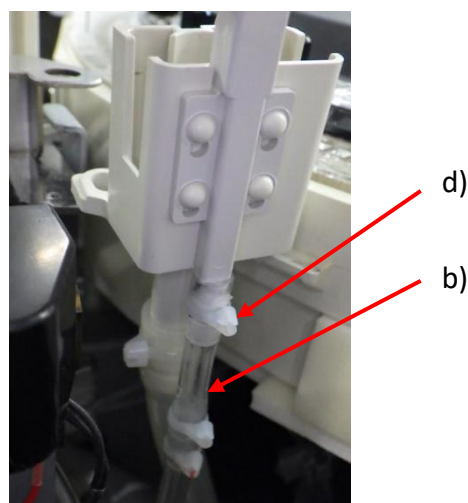


Fig.20 Supply tube

a) CLEANING BATH S1 A (replacement part)

b) Supply tube

c) Drain tube

d) Wire binder

(13) Fix the CLEANING BATH S1 A and SPACER with the screws (M644039). (2 screws)

The screws were tightened temporarily.

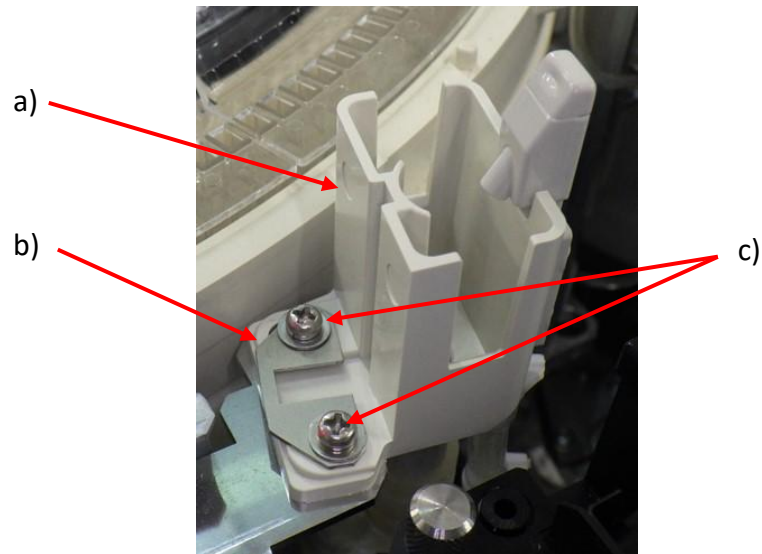


Fig.21 CLEANING BATH S1 A

a) CLEANING BATH S1 A

b) SPACER

c) Screw(M644039)

(14) Insert the SPACER10 between CLEANING BATH S1 A and base of Rinse Station.

Be careful not to drop the spacer and insert it all the way in. Tighten the screws to secure the CLEANING BATH S1 A.

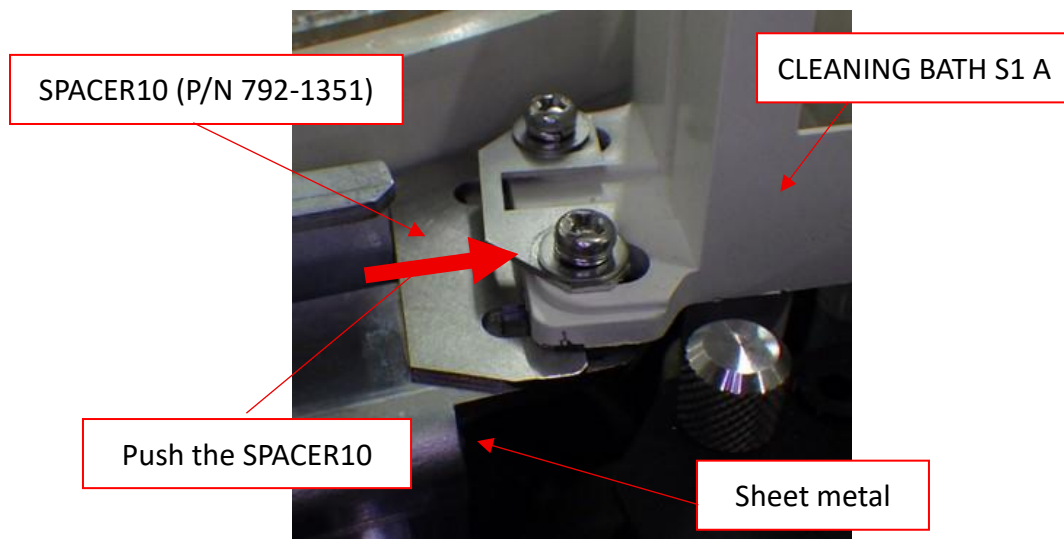


Fig.22 SPACER10

(15) Attached the WS FIXING PLATE 3 (P/N792-2777). (1 screw)

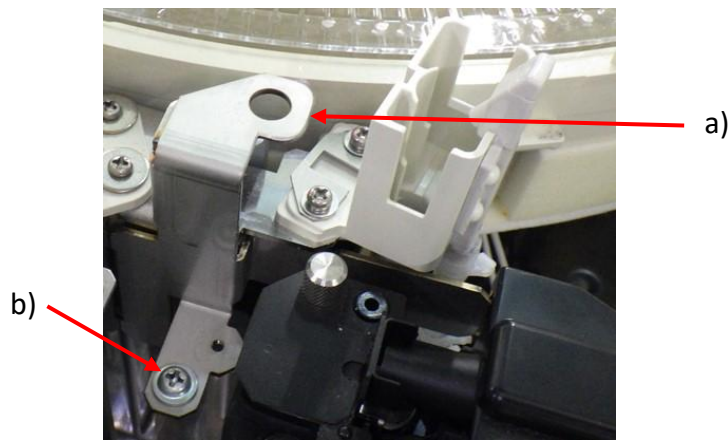


Fig.23 WS FIXING PLATE 3

a) WS FIXING PLATE 3

b) Screw

(16) Attached the PIPETTER S1 ASSY removed in step (8).

(17) Reconnect the connector and nipple that removed in step (8).

Make sure the S1 probe doesn't hit the CLEANING BATH S1 A when the S1 arm rotates by hand. If the S1 probe hit the CLEANING BATH S1 A, shift the position of the CLEANING BATH S1 A.

(18) Remove 3 covers.

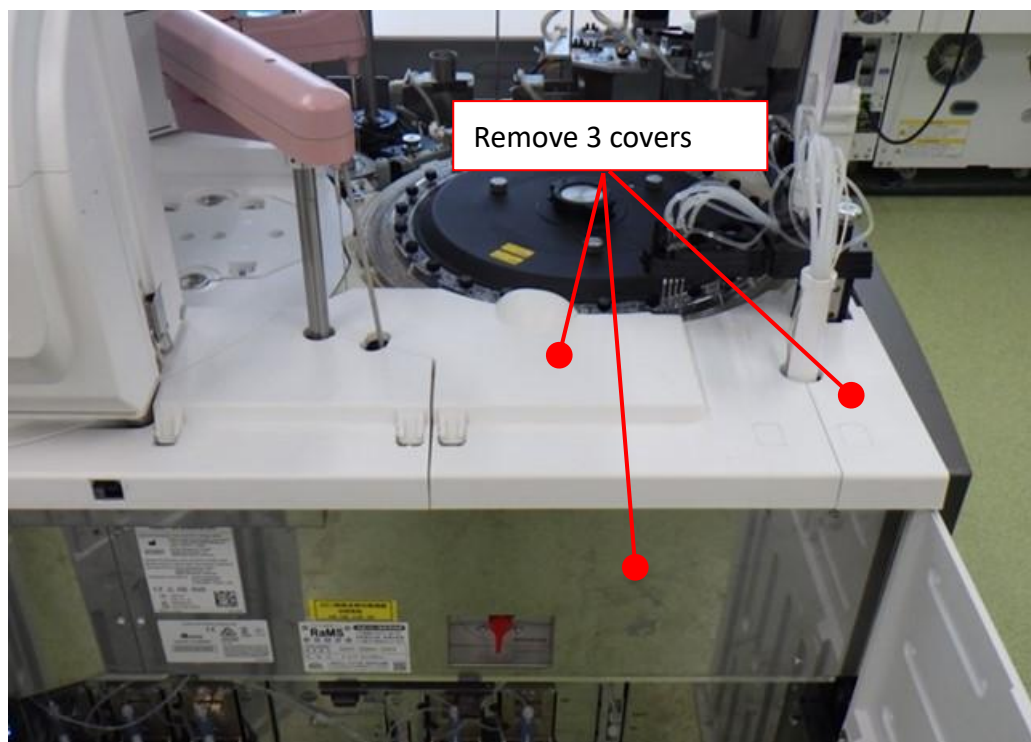


Fig. 24 Covers

(19) Remove a screw for fixing S1 needle valve.

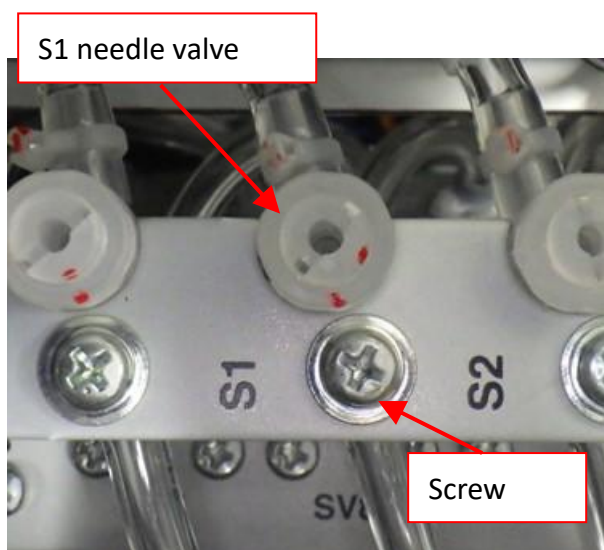


Fig.25 S1 needle valve

- (20) Cut off Wire Binder and remove a Tube. When removing the tube, water may spill out from the tube, put a paper towel under the tube.

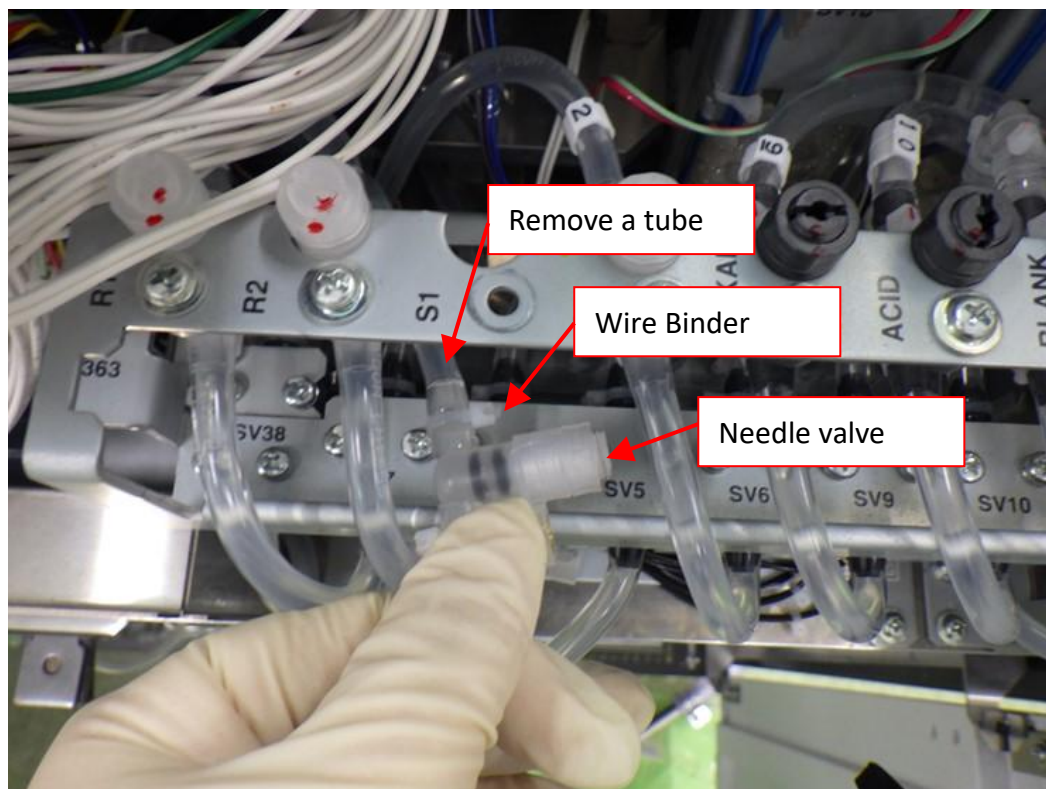


Fig. 26 Remove S1 needle valve

- (21) Connect the TUBE NV2S1RS ASSY (P/N 792-7754). If the tip of the tube is deformed, cut it off around 10mm from tip. Secure each end of the TUBE NV2S1RS ASSY with 2 Wire Binds (J839050).

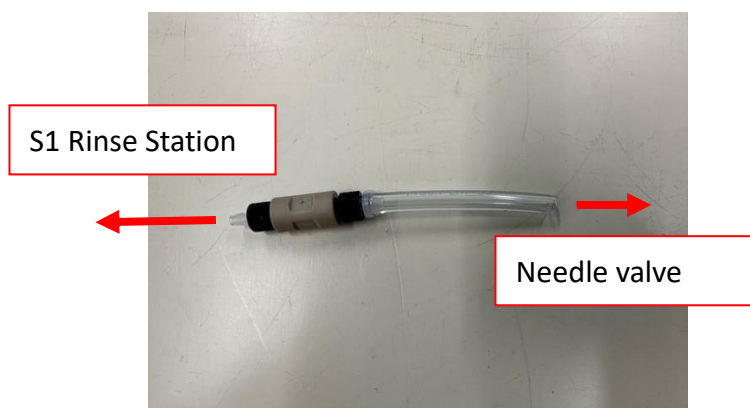


Fig. 27 TUBE NV2S1RS ASSY

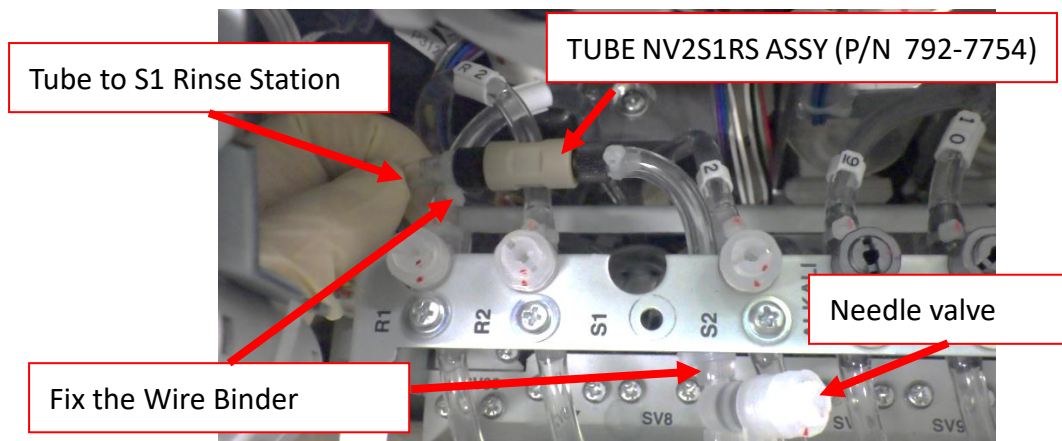


Fig. 28 Connect TUBE NV2S1RS ASSY

(22) Remove the Screw.

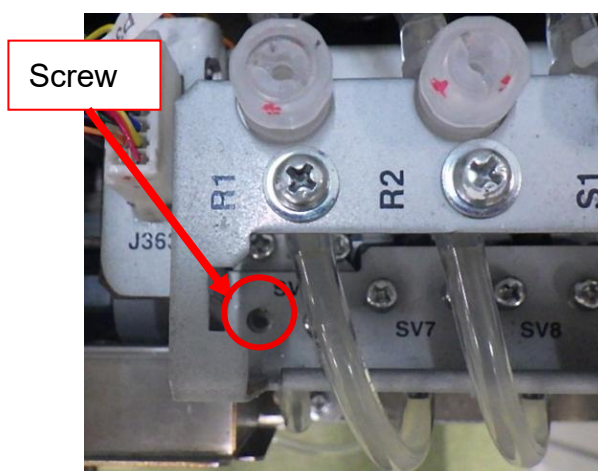


Fig.29 Remove Screw

(23) Installing the JOINT FIX PLATE ASSY (P/N 792-7756) with M641135.

JOINT FIX PLATE ASSY is installed so that a part of it is under sheet metal of needle valve.

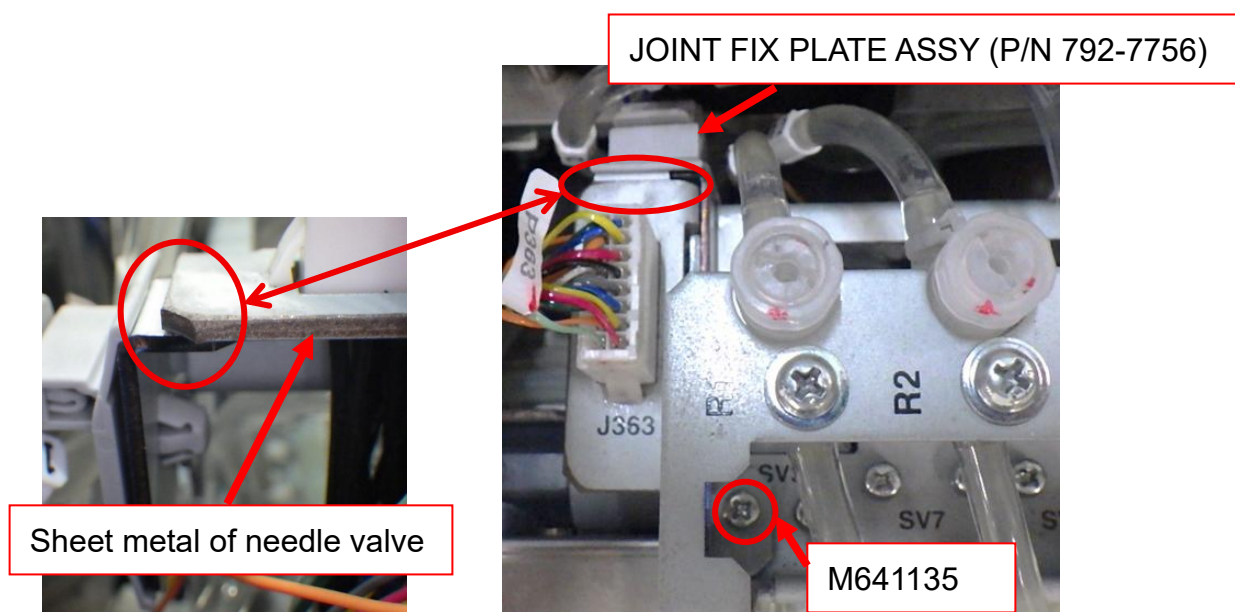


Fig. 30 JOINT FIX PLATE ASSY

(24) Return the needle valve to sheet metal and fixed by screw.

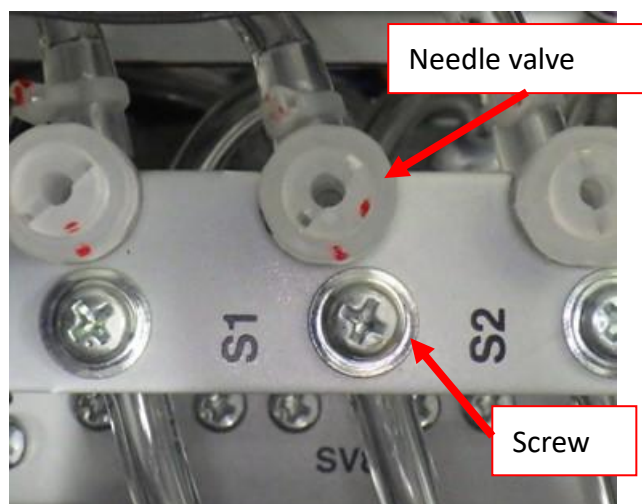


Fig. 31 Return the needle valve

- (25) Fix the TUBE NV2S1RS ASSY (P/N 792-7754) to the clamp of the JOINT FIX PLATE ASSY (P/N 792-7756). Make sure the cable is not pushing against the tube and that the tube does not kinking.

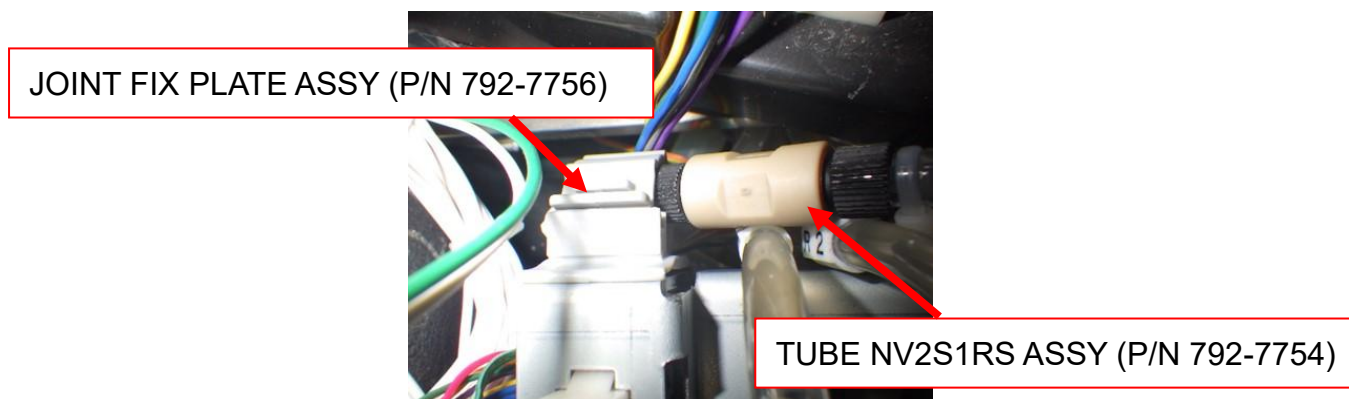


Fig.32 Clamping

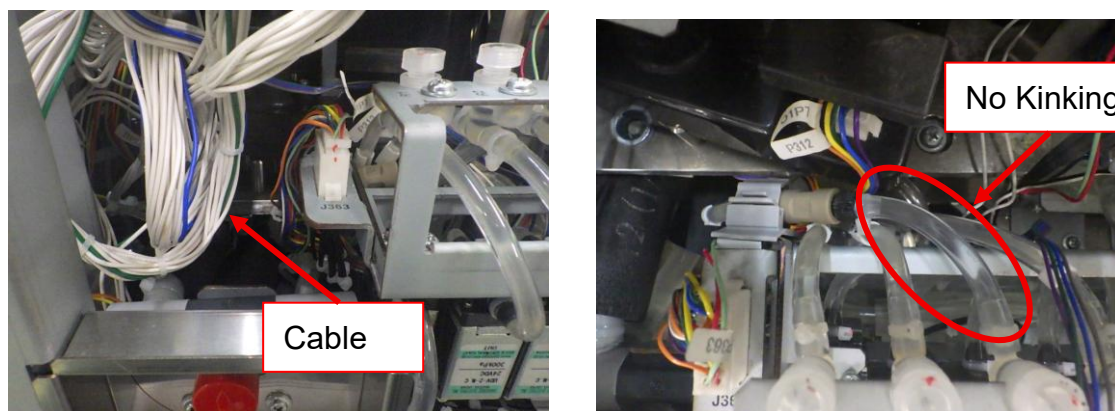


Fig. 33 Tube Routing

- (26) Replacement parts from the TB cover 5 to the COVER TB5A (P/N 792-2881).
 (27) Power on the instruments.

Adjust the horizontal (Washing position) of CLEANING BATH S1 A and Sample prob S1.

Adjust the outside rinse water volume (Sample probe S1). See below for adjustment procedures.

Instrument components / Fluidics / Water volume / Check / adjust / Adjusting the cell rinse water volume

(28) Attached Covers and Parts removed in step (3)-(6), (17).

Install COVER TB5A (P/N 792-2881) and COVER TB7 ASSY (P/N 792-0860).

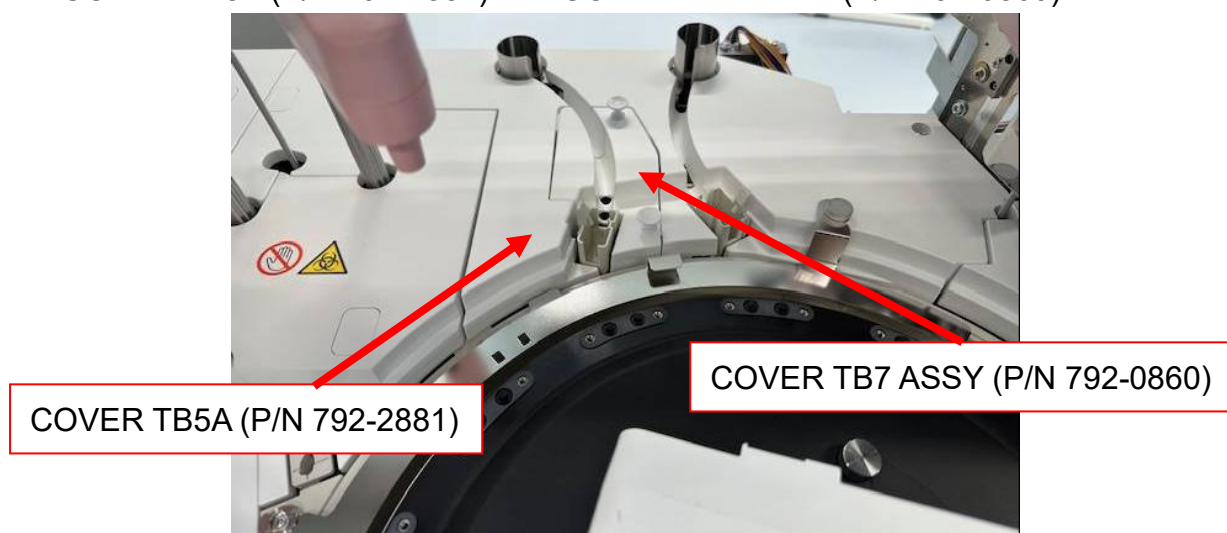


Fig. 34 Install COVER

end of document