


cobas[®] pure – c303 Adjustments



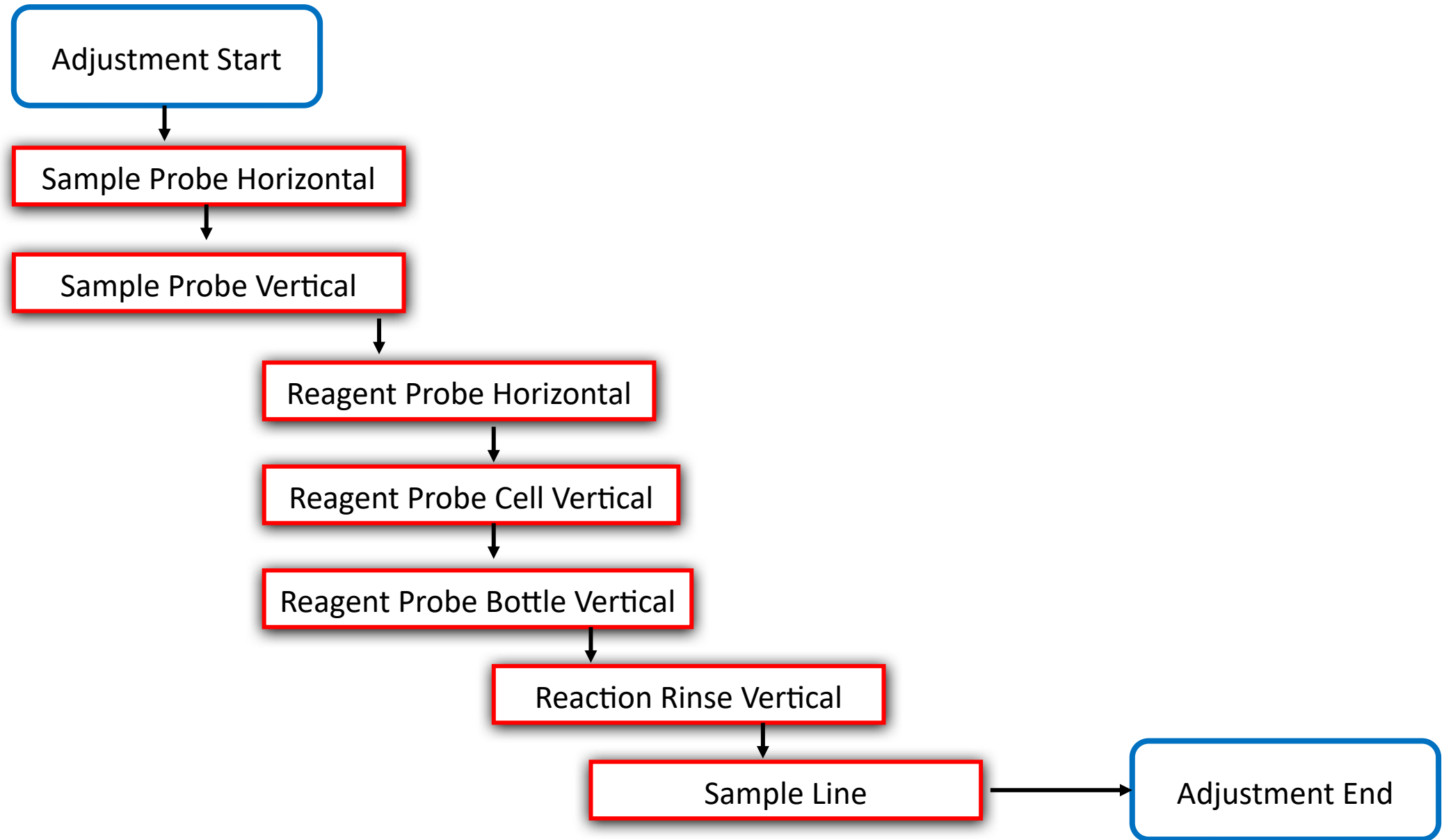

Sample Probe Horizontal

- Preparation :

Tool 1 (S Nozzle Jig) is needed.

P/N	Tool (1)	
	798-3280	798-3281
Item name	S Nozzle Jig Assy	S Nozzle Jig Assy
Appearance		

Remove the cover shown in Figure 10 in advance.


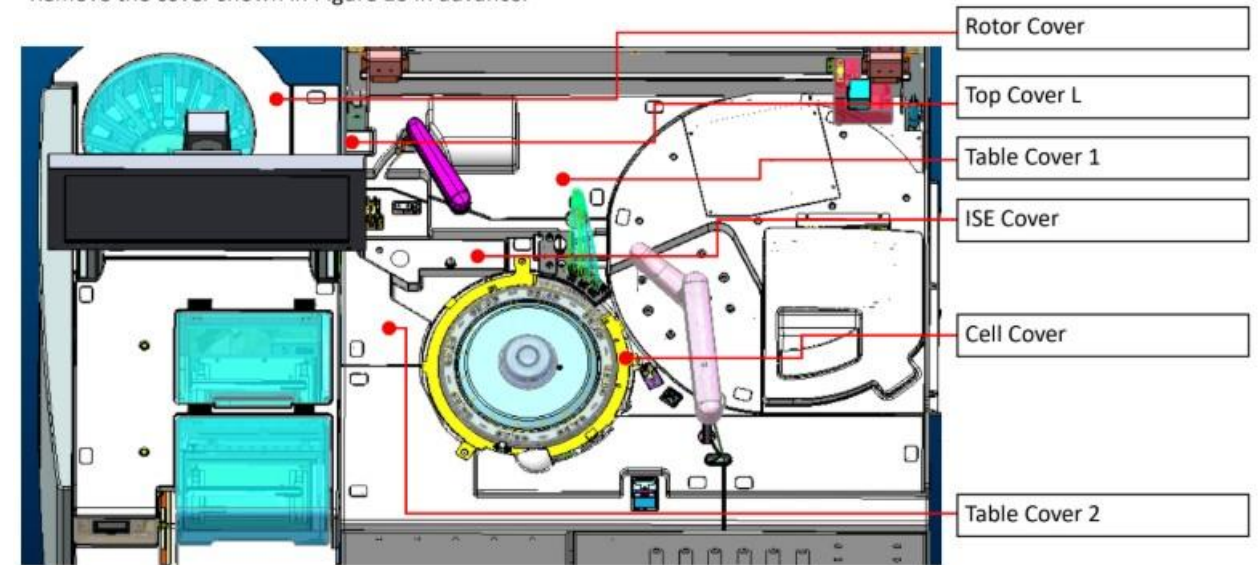


Figure 10. Removing the Cover


Washing (1/7)




1. Adjust the Probe , so that it is at the center of washing station.
2. Press 'Accept & Next' after adjustment.

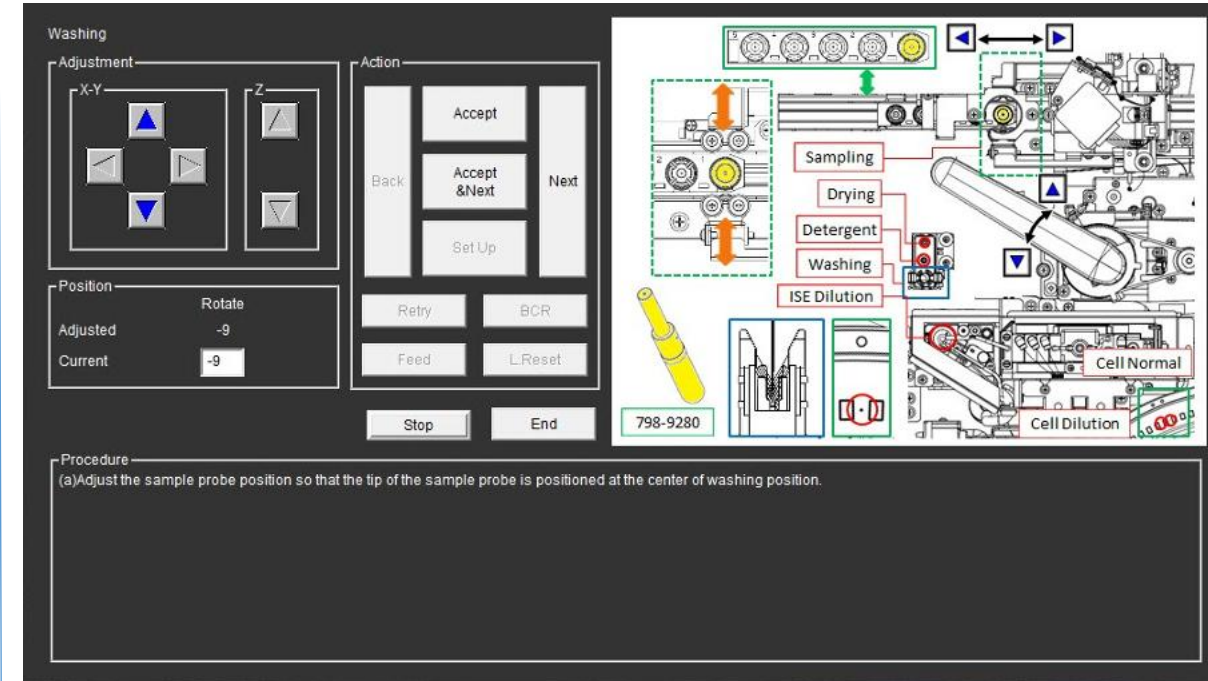

Vacuum (Drying) (2/7)

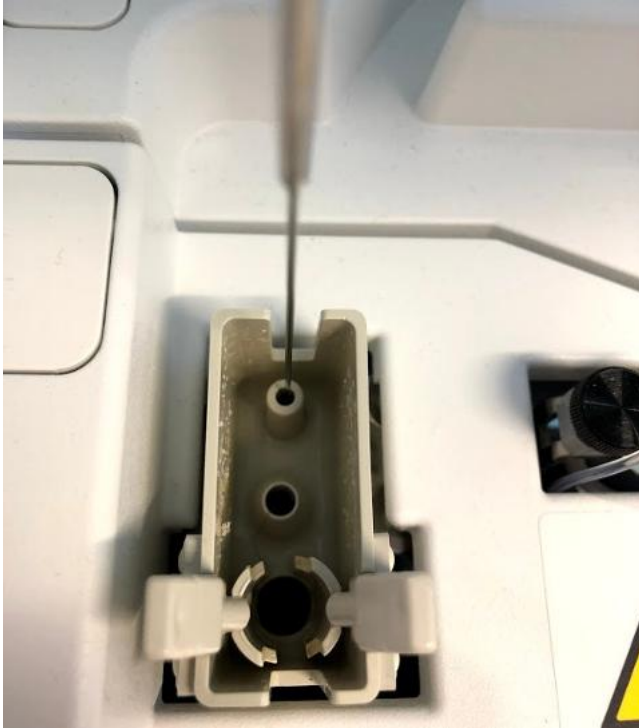

1. Adjust the Probe , so that it is at the center of vacuum hole.
2. Press 'Accept & Next' after adjustment.

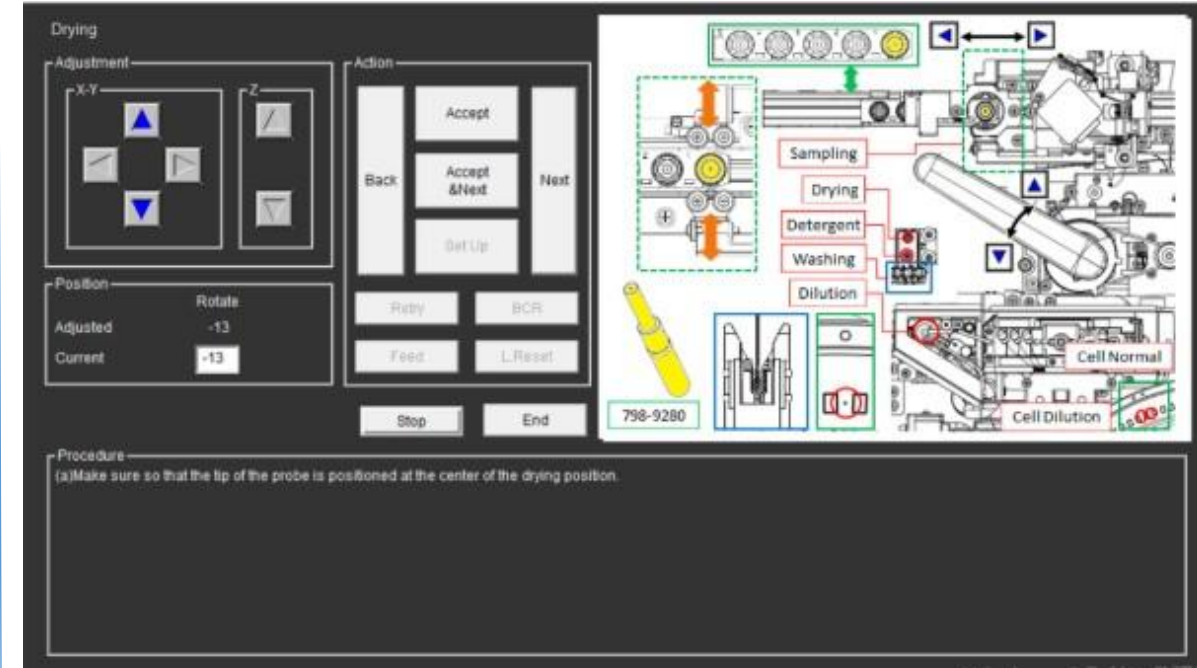

Sampling (3/7) – Preparation :



Place the adjustment tool on pos.1 on a Sample Rack and set the rack at Sample Line of c303. Press 'Setup' to move the rack to Sampling Position.

Sampling

Adjustment

X-Y: [Up Arrow] [Down Arrow] [Left Arrow] [Right Arrow]

Z: [Up Arrow] [Down Arrow]

Position

S Line: Adjusted 3, Current 3

Rotate: Adjusted -33, Current -33

Action

Back, Accept, Accept & Next, Set Up, Next, Retry, BCR, Feed, L.Reset, Stop, End

Procedure

- Place an adjustment tool (P/N: 798-9280) on the rack Pos. 1. Then place the rack with the adjustment tool on the sample line.
- Press the "Feed" button to move the rack to SP.
- Make sure so that the tip of the probe fits into the hole of the adjustment tool at SP.
- Adjust the active gripper position so that the hole of the adjustment tool is positioned at the same place during open / close the active gripper. Press the "Set Up" button to open / close the active gripper.
- After making sure, press the "Feed" button repeatedly until the rack moves to the setting position.
- Take out the rack from the setting position.

798-9280

Sampling Process Diagram: Sampling, Drying, Detergent, Washing, Dilution, Cell Normal, Cell Dilution

Sampling (3/7)



1. Adjust the Probe , so that it is at the center of the adjustment tool.
2. Press 'Setup' to Open/Close active Gripper and check whether the position of the adjustment changes.
3. Press 'Accept' after adjustment.
4. Press 'Feed' several times to get the rack out to setting position. Remove rack.
5. Press 'Next'.

Sampling

Adjustment

X-Y: [Up Arrow] [Down Arrow] [Left Arrow] [Right Arrow]

Z: [Up Arrow] [Down Arrow]

Position

S Line: Adjusted 3, Current 3

Rotate: Adjusted -33, Current -33


Action

Buttons: Back, Accept, Accept & Next, Next, Set Up, Retry, BCR, Feed, L.Reset, Stop, End


Procedure

(a) Place an adjustment tool (P/N: 798-9280) on the rack Pos. 1. Then place the rack with the adjustment tool on the sample line.
(b) Press the "Feed" button to move the rack to SP.
(c) Make sure so that the tip of the probe fits into the hole of the adjustment tool at SP.
(d) Adjust the active gripper position so that the hole of the adjustment tool is positioned at the same place during open / close the active gripper. Press the "Set Up" button to open / close the active gripper.
(e) After making sure, press the "Feed" button repeatedly until the rack moves to the setting position.
(f) Take out the rack from the setting position.


Cell Normal (4/7)

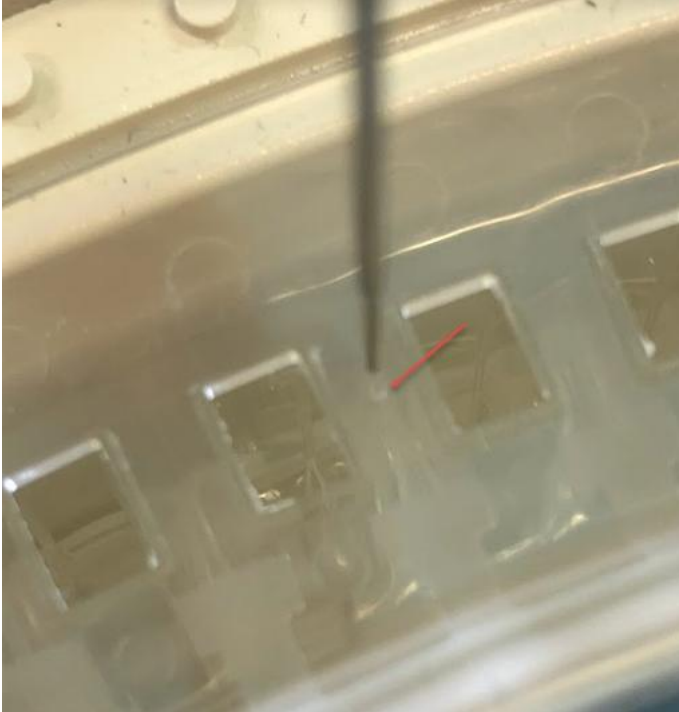


1. Adjust the Probe , so that it is positioned at the center of the spot as shown.
2. Press 'Accept & Next' after adjustment.

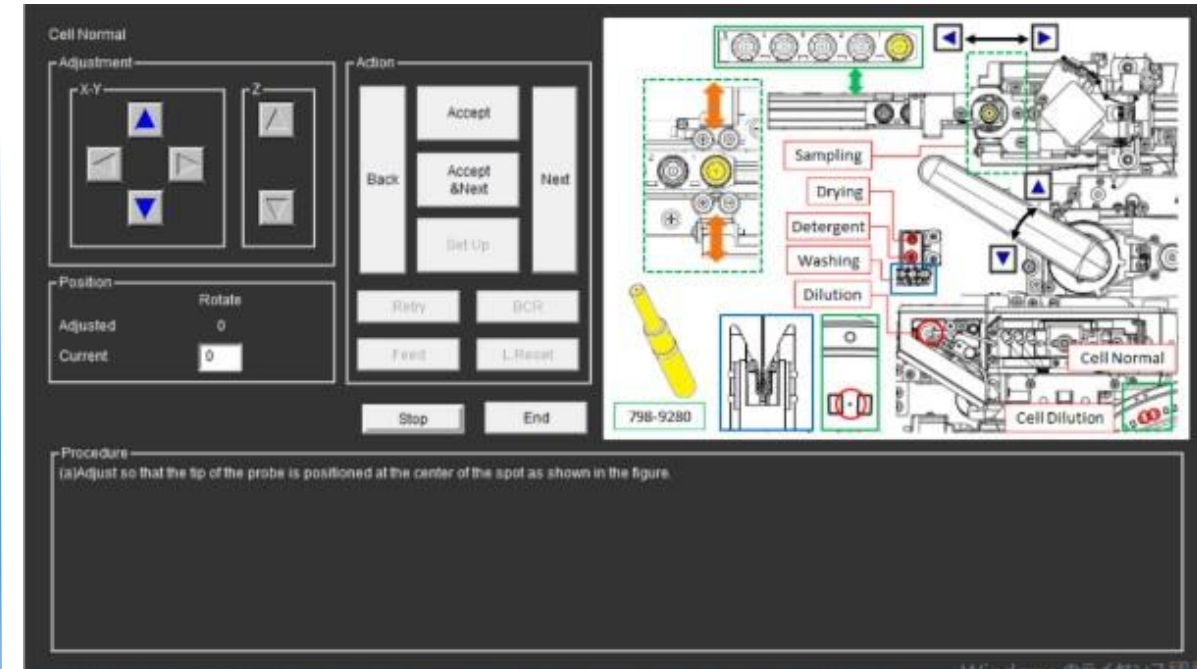

Cell Dilution (5/7)

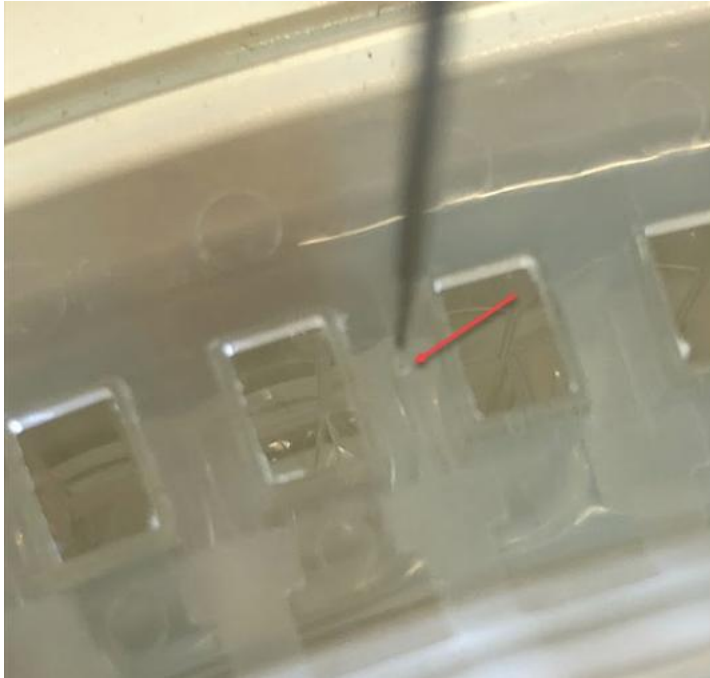


1. Adjust the Probe , so that it is positioned at the center of the spot as shown.
2. Press 'Accept & Next' after adjustment

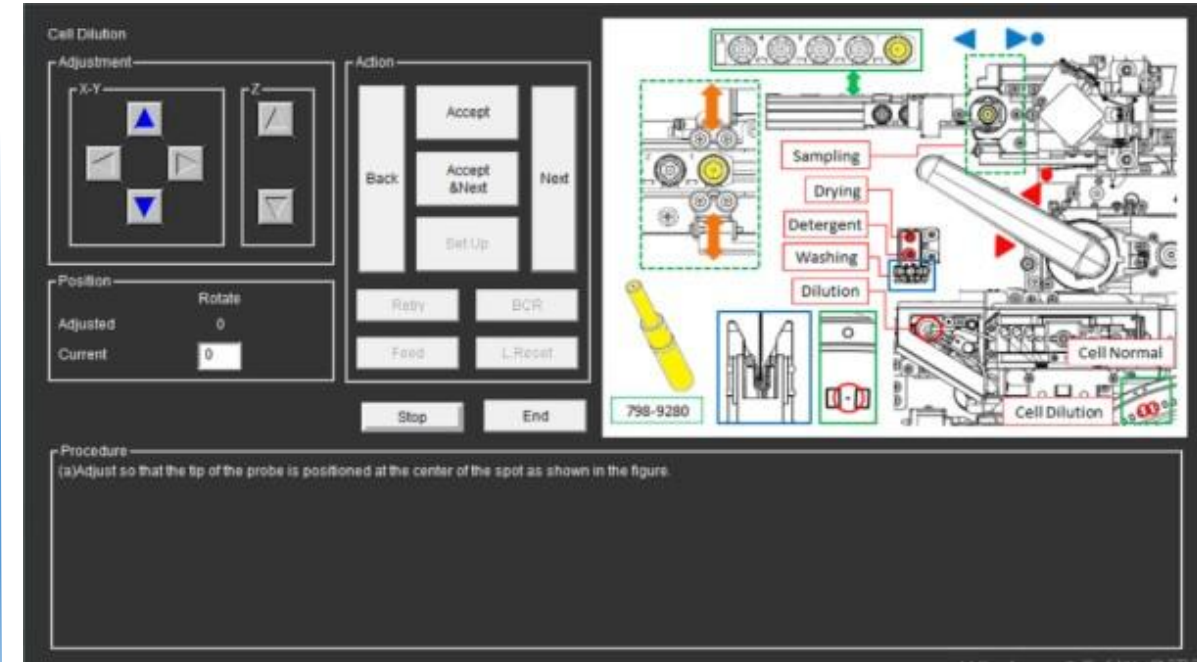

ISE Dilution Vessel (6/7)




1. Adjust the Probe , so that it is positioned at the center of the dilution vessel as shown.
2. Press 'Accept & Next' after adjustment.

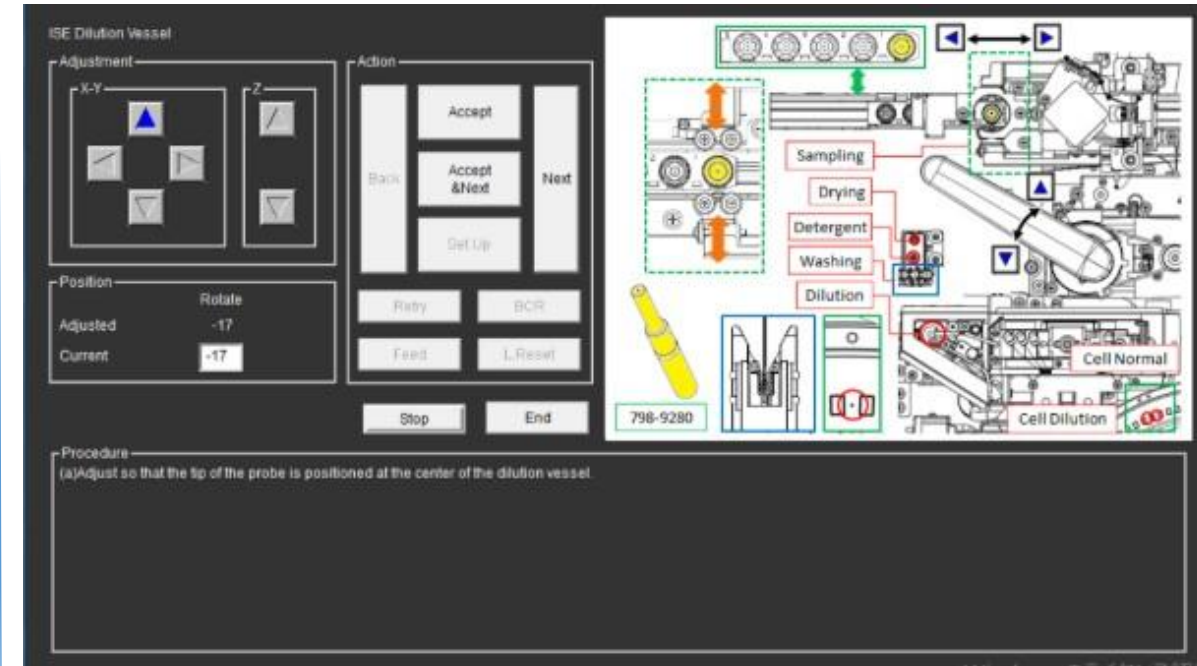

SWS (Detergent)(7/7)




1. Adjust the Probe , so that it is positioned at the center of the detergent hole as shown.
2. Press 'Accept & Next' after adjustment.

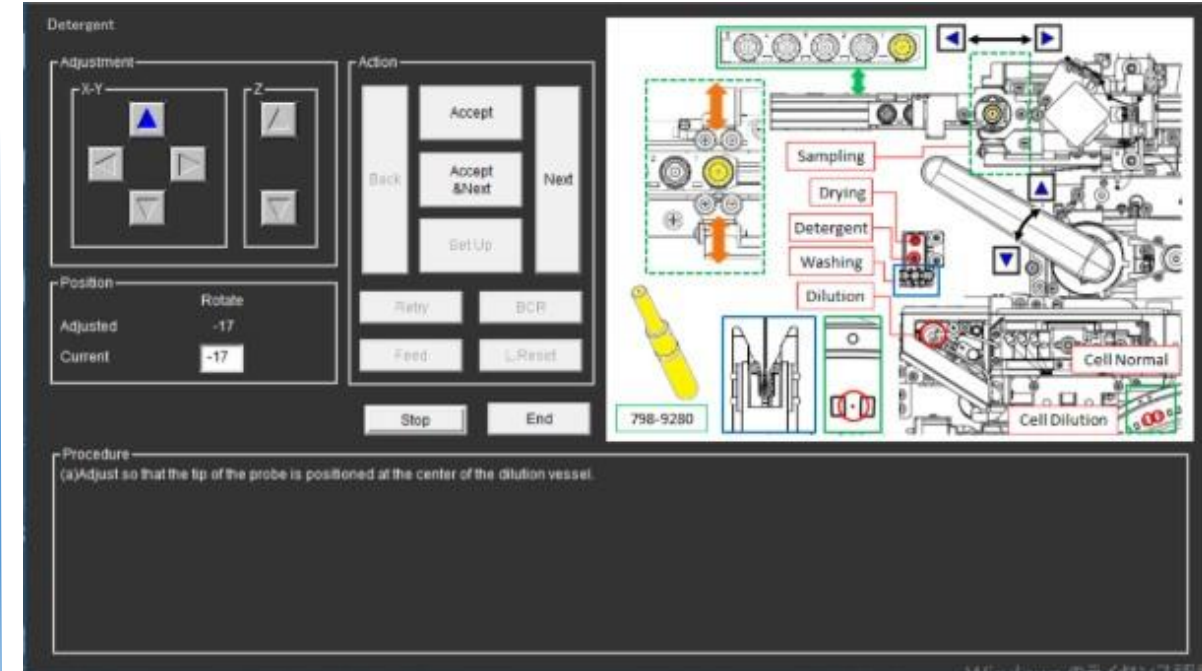

Sample Probe Vertical

- Preparation :
Tool 2 (S.P HM2 Jig) is needed.


	Tool (2)
P/N	792-6427
item name	S.P Jig HM2
Appearance	

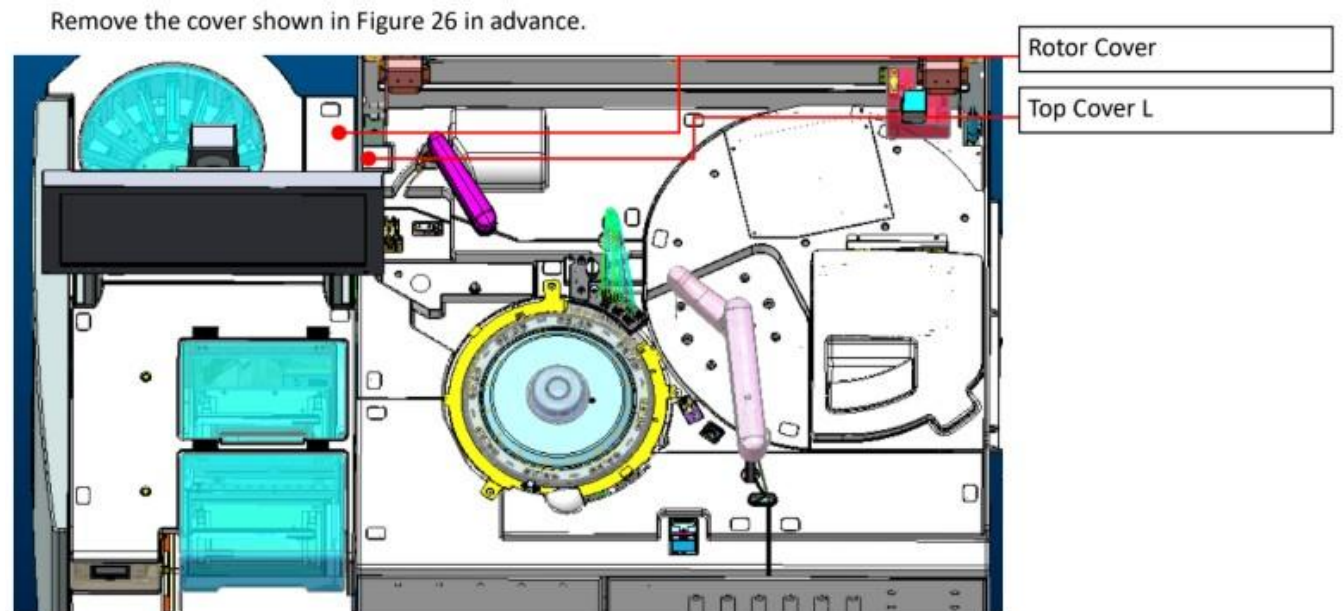

Sample Probe Vertical (1/1) – Preparation :



Place the height adjustment tool on pos.1 on a Sample Rack and set the rack at the Sample Line of c303.

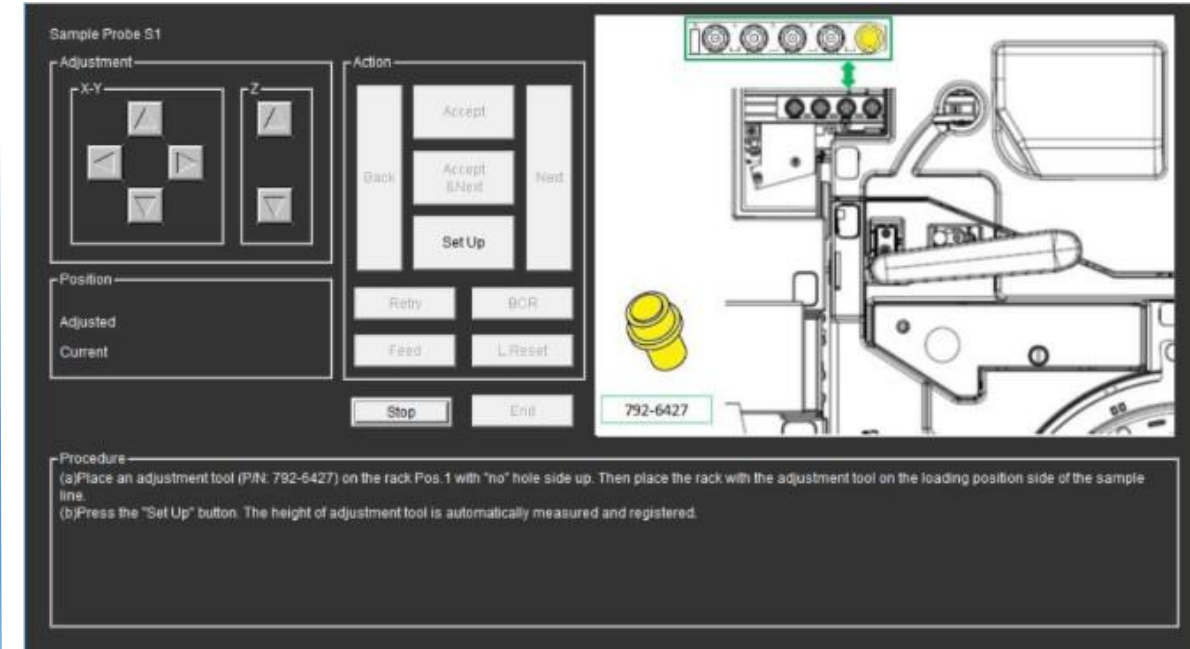


Adjustment Flowchart – Start Page

Sample Probe Vertical (1/1)



1. Press 'Setup' .
2. Probe will move automatically , first to Sampling Position , then to cell bottom and finally ISE dilution vessel bottom ..
3. Press 'End' to finish.

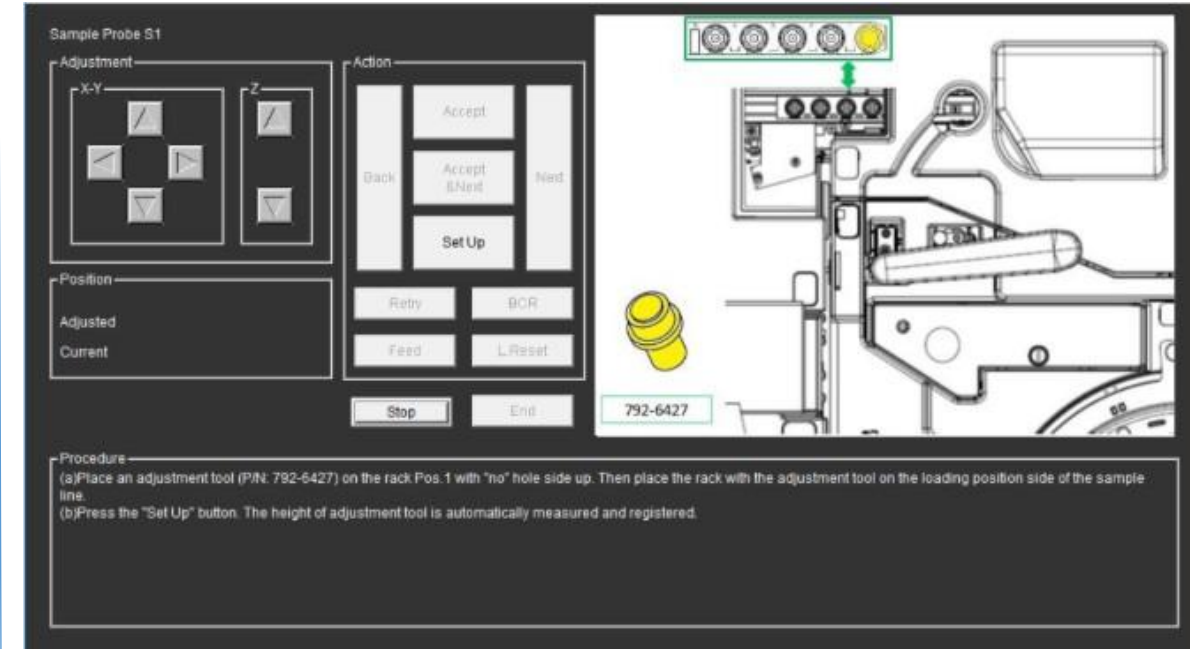

Reagent Probe Horizontal

- Preparation :

Remove the cover shown in Figure 30 in advance.



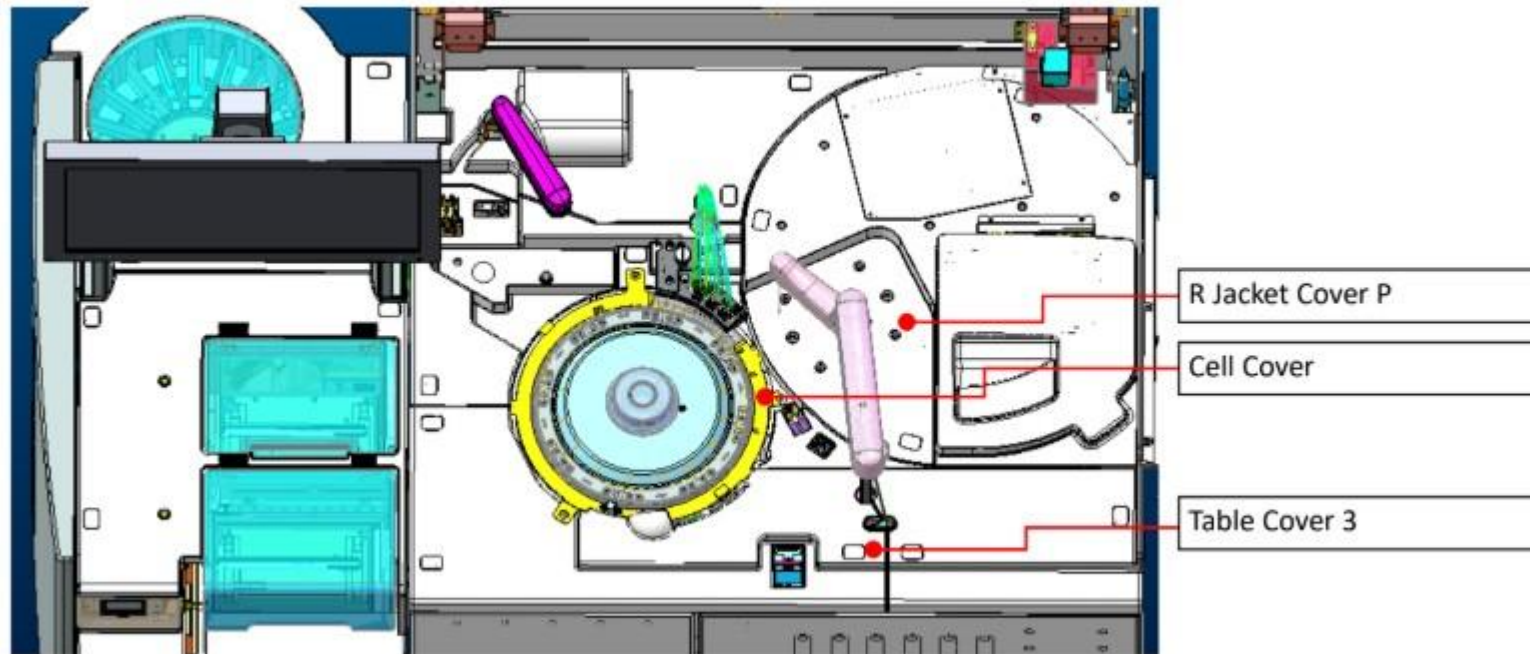


Figure 30. Removing the Cover


Washing (1/14)

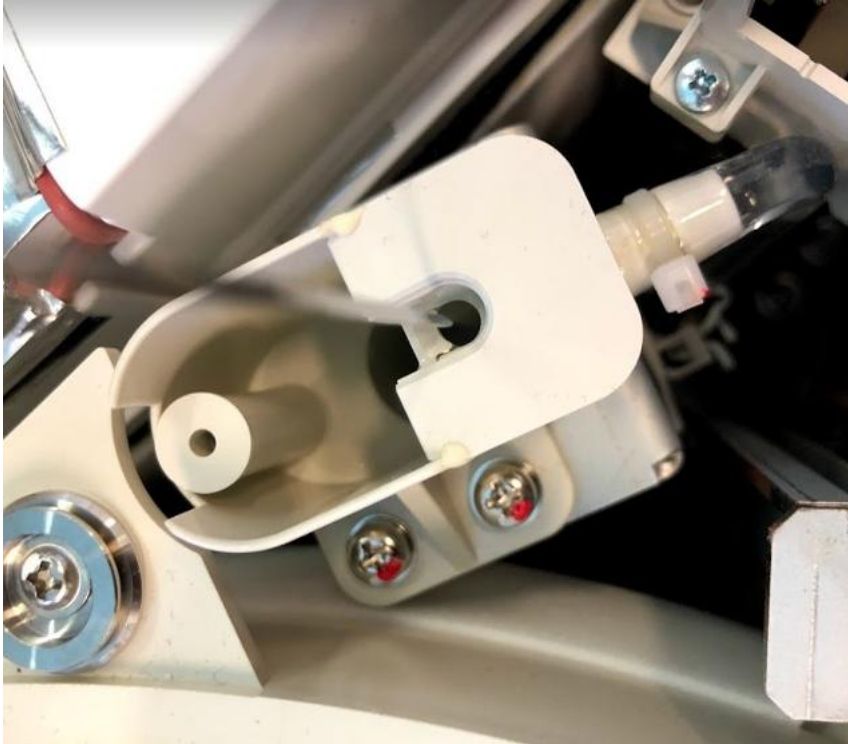


1. Adjust the Probe , so that it is at the center of Washing hole.
Move up&down to check the position.
2. Press 'Accept & Next' after adjustment.

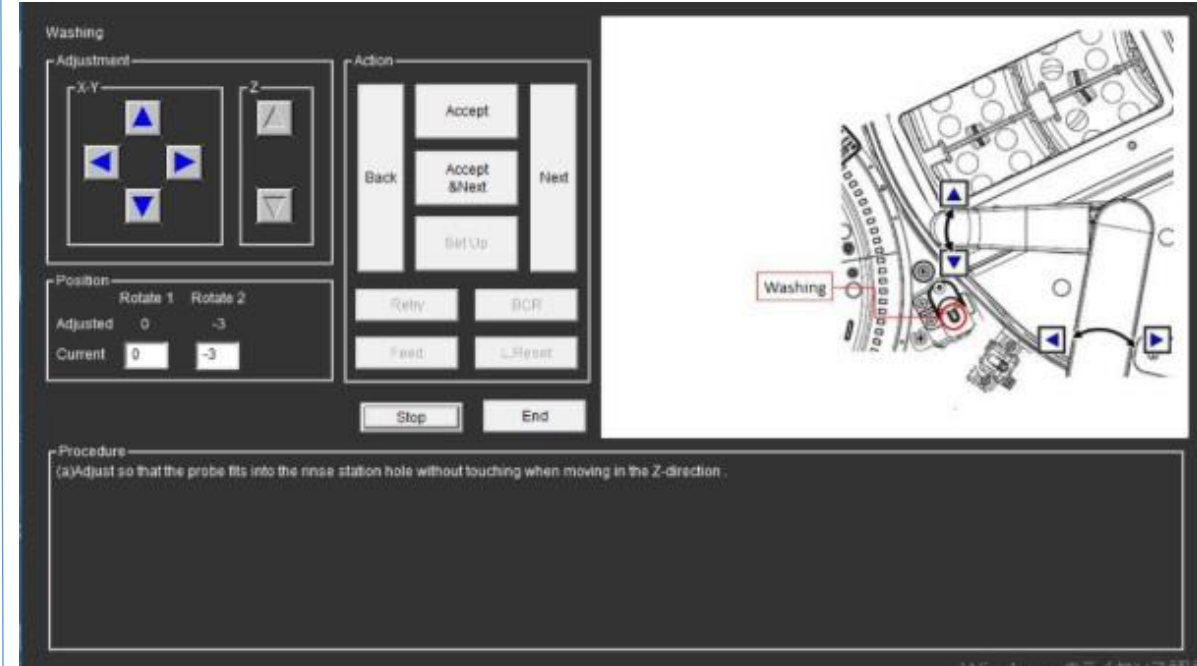


Drying (2/14)

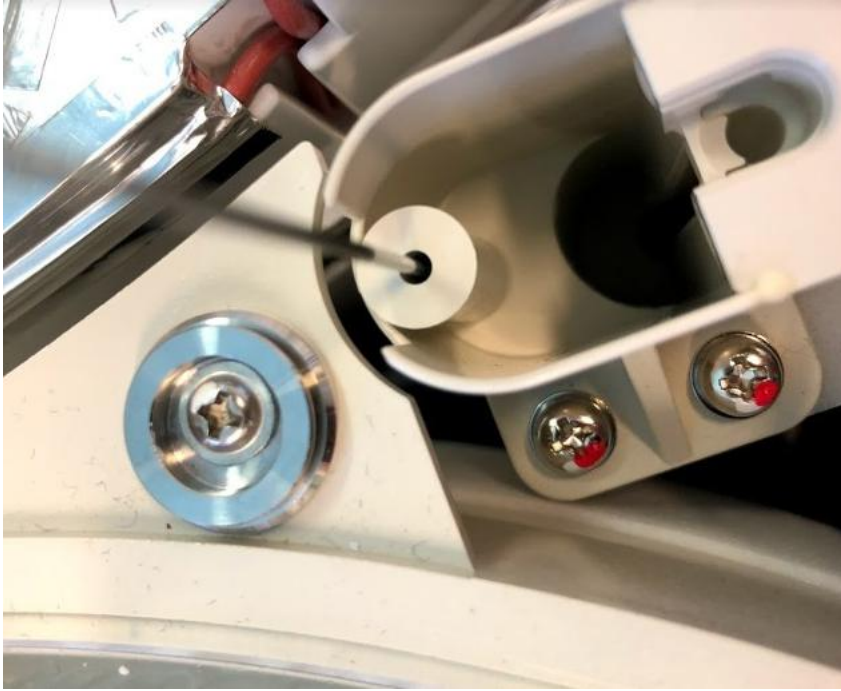


1. Adjust the Probe , so that it is at the center of Drying hole. Move up&down to check the position.
2. Press 'Accept & Next' after adjustment.

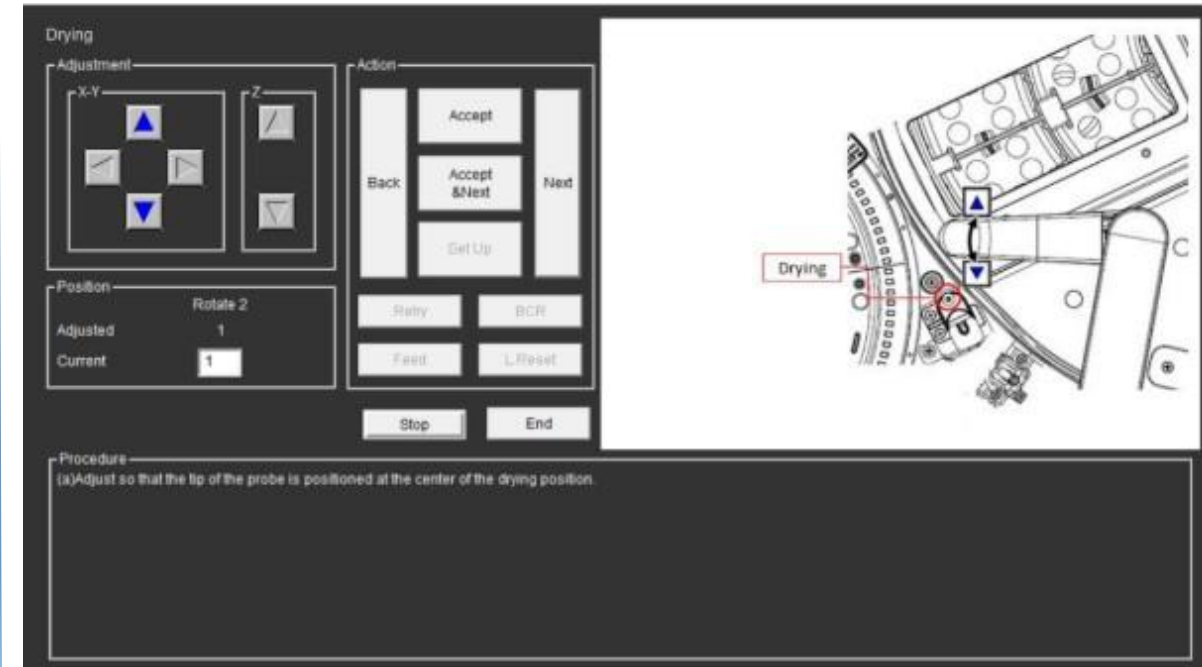

R Disk CCW (3/14)

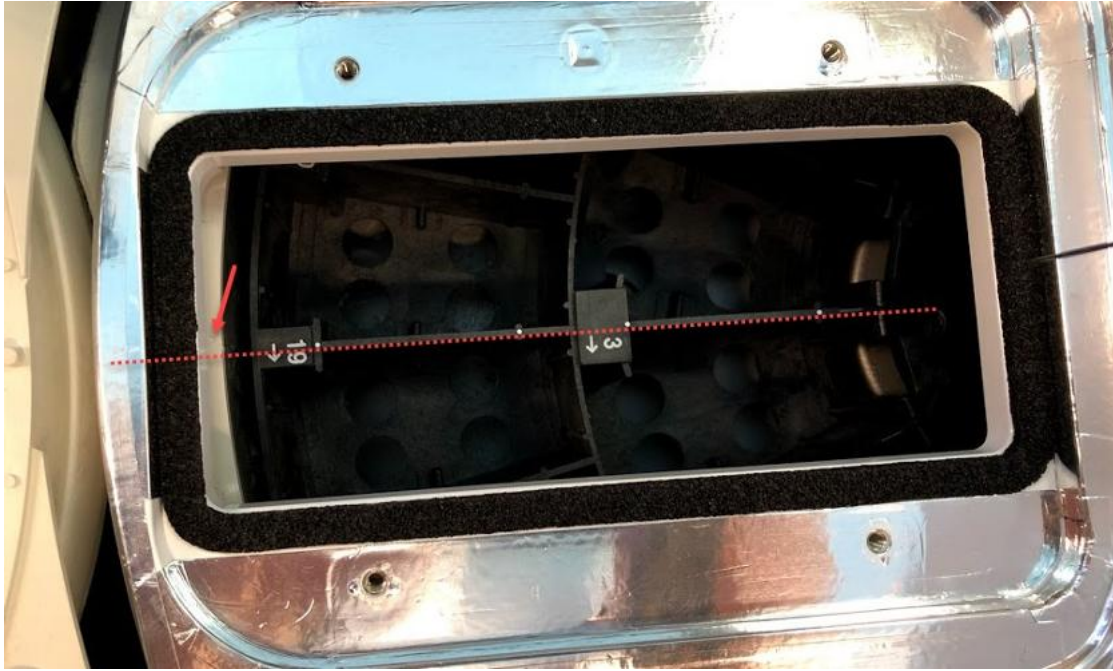


1. Adjust the Reagent Disk position , so that the indicator on the disk cover and the inner line are aligned .
2. Press 'Accept & Next' after adjustment.

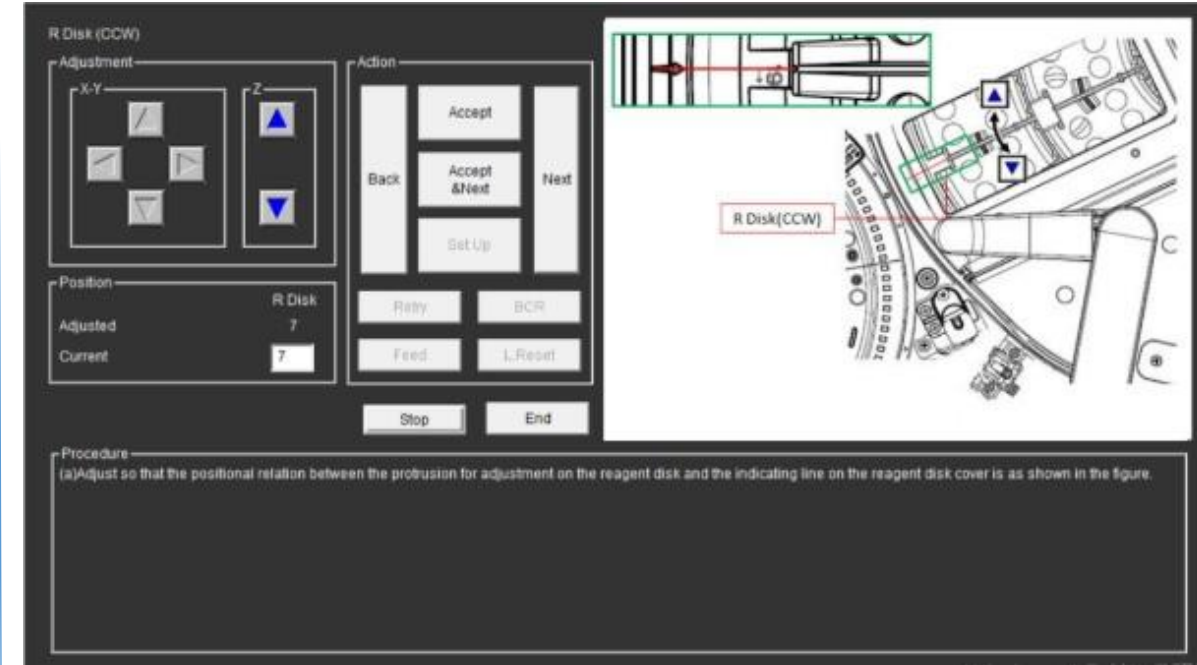


Position 1 (4/14)




1. Adjust the probe , so that the tip of the probe is centered on the position.
2. Press 'Accept & Next' after adjustment.

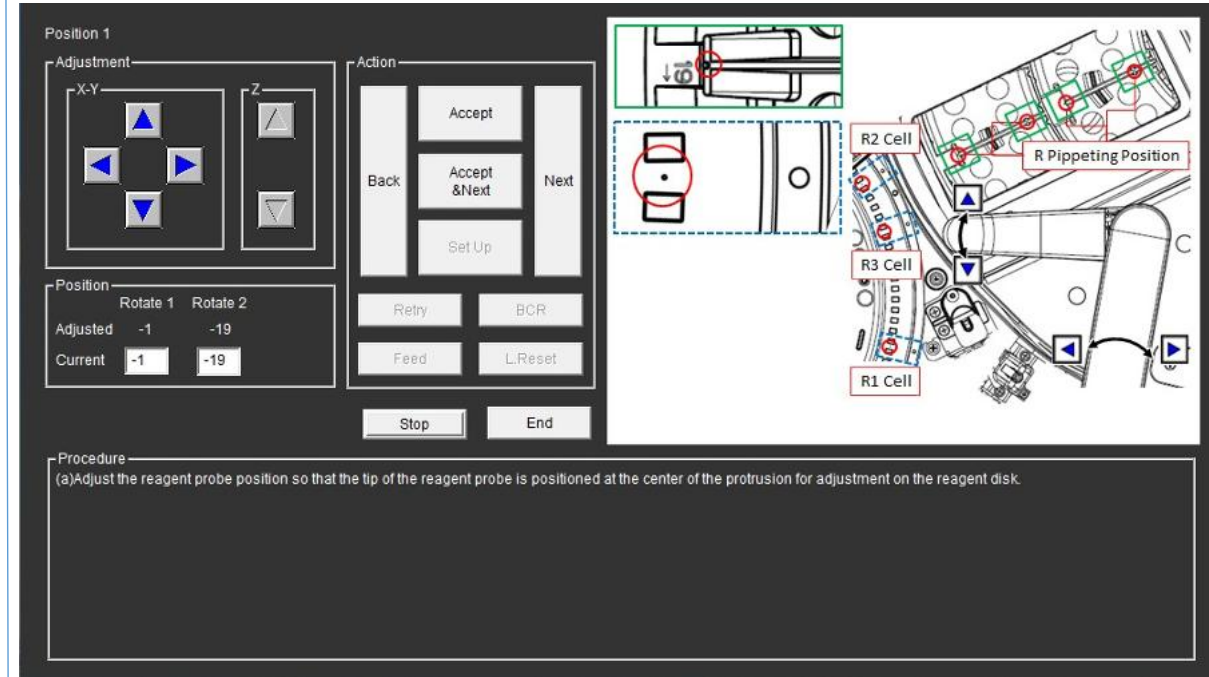

R1 Cell (5/14)

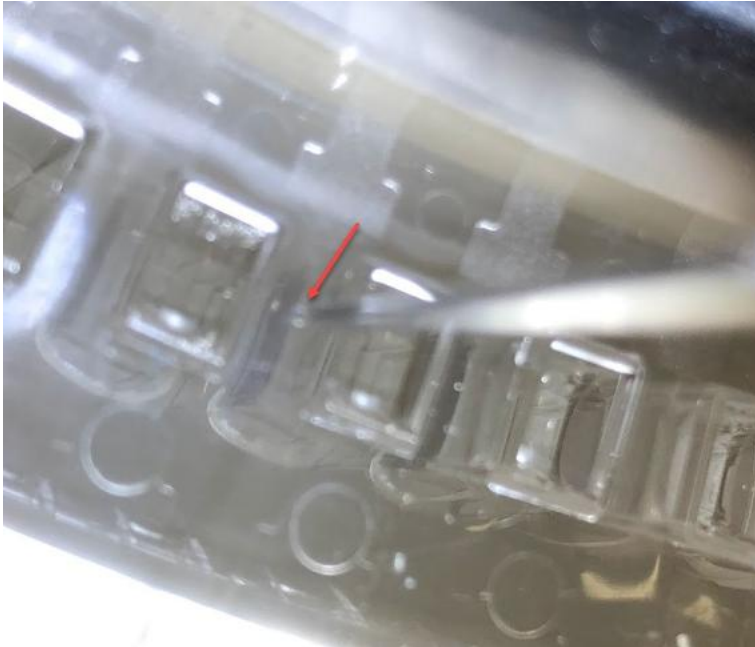


1. Adjust the Probe , so that it is positioned at the center of the spot as shown.
2. Press 'Accept & Next' after adjustment

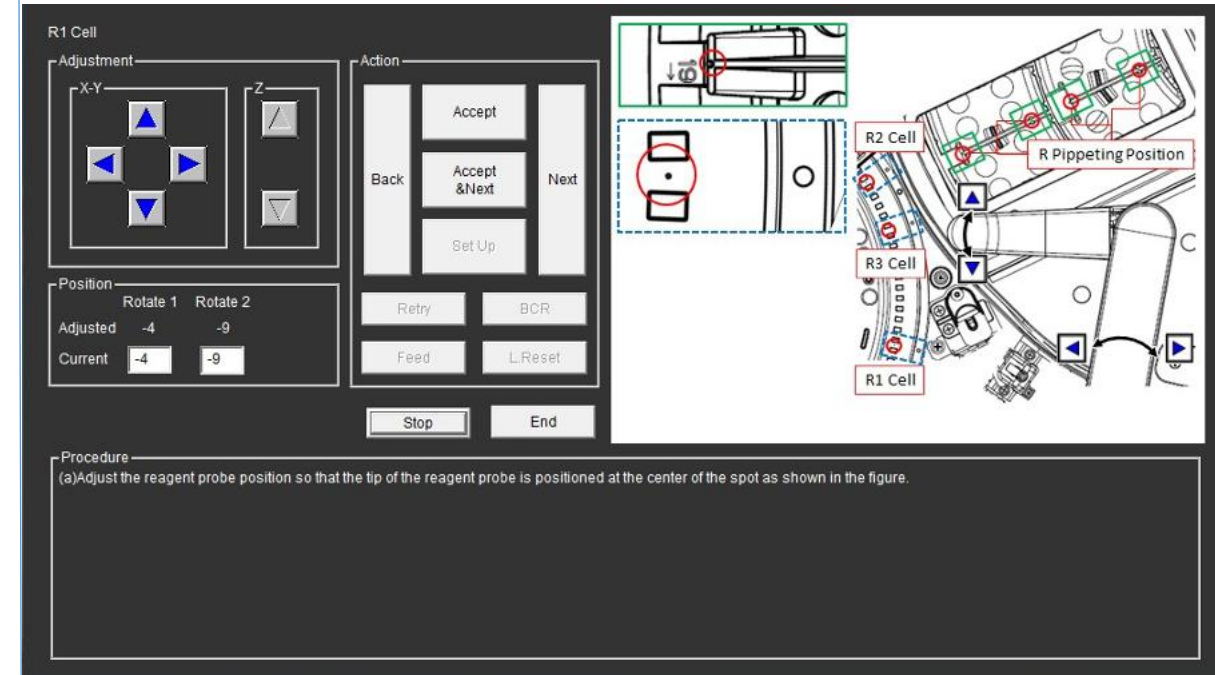

Position 3 (6/14)




1. Adjust the probe , so that the tip of the probe is centered on the position.
2. Press 'Accept & Next' after adjustment.

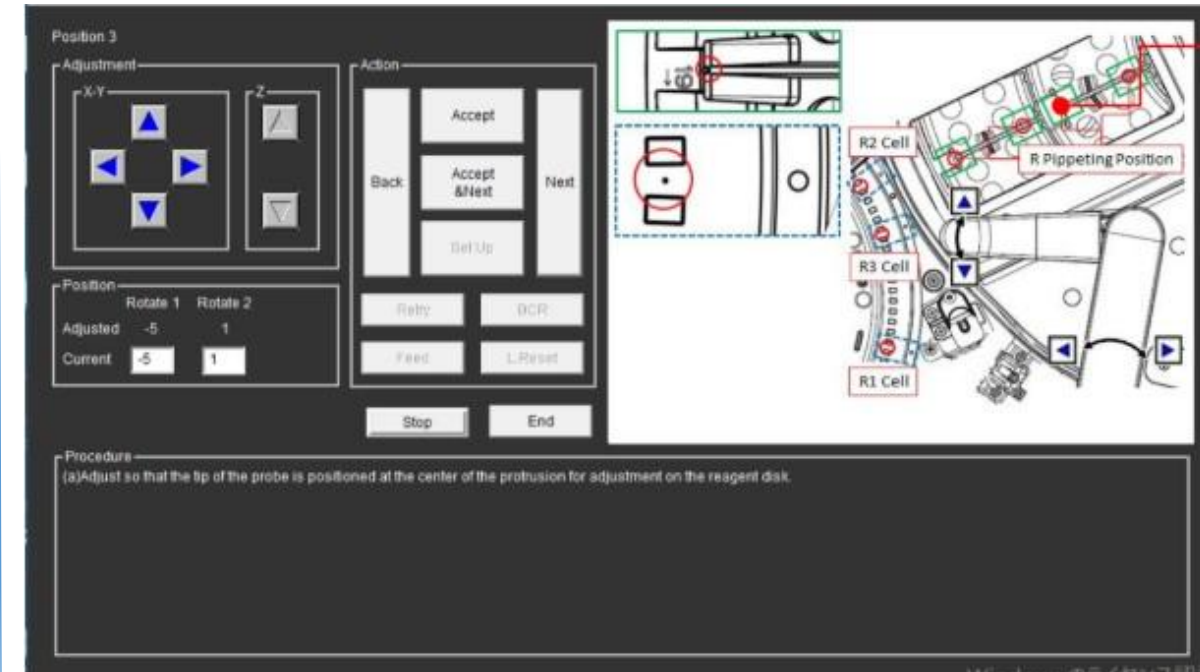

R2 Cell (7/14)




1. Adjust the Probe , so that it is positioned at the center of the spot as shown.
2. Press 'Accept & Next' after adjustment

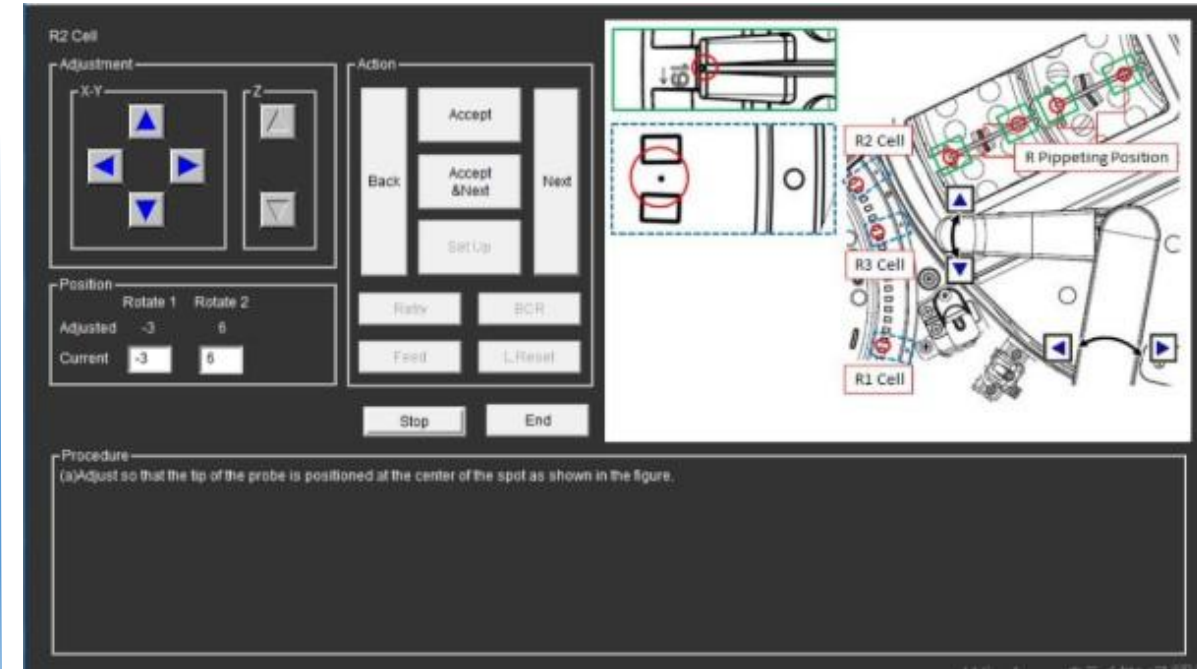

Position 4 (8/14)




1. Adjust the probe , so that the tip of the probe is centered on the position.
2. Press 'Accept & Next' after adjustment.

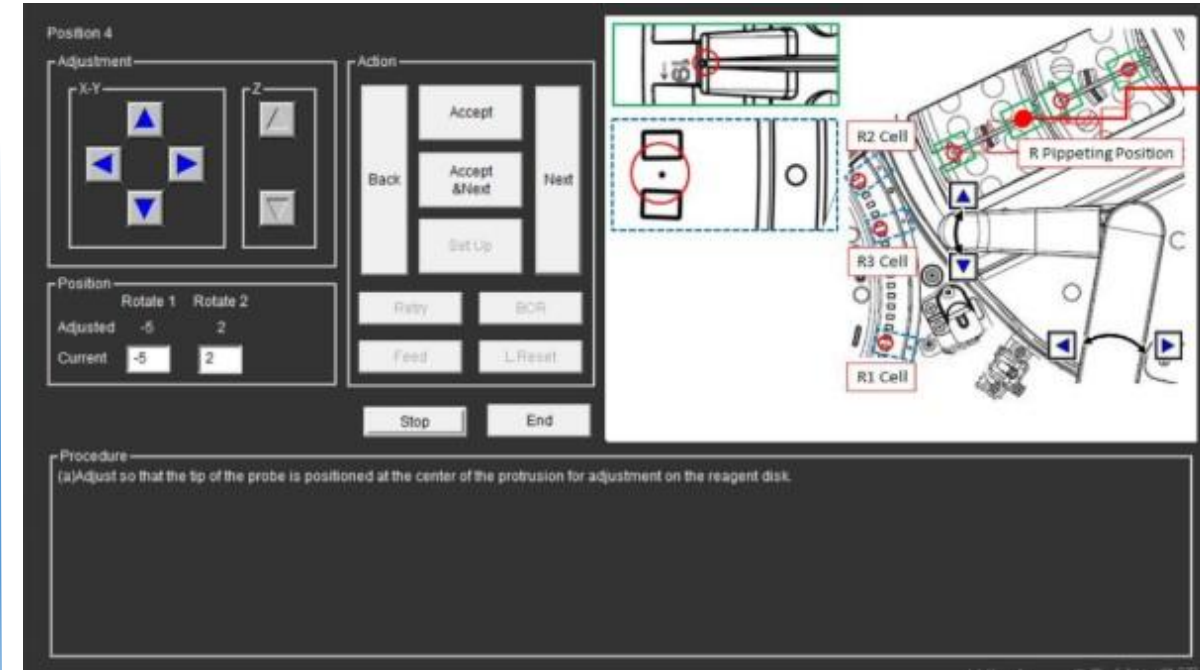


R3 Cell (9/14)




1. Adjust the Probe , so that it is positioned at the center of the spot as shown.
2. Press 'Accept & Next' after adjustment

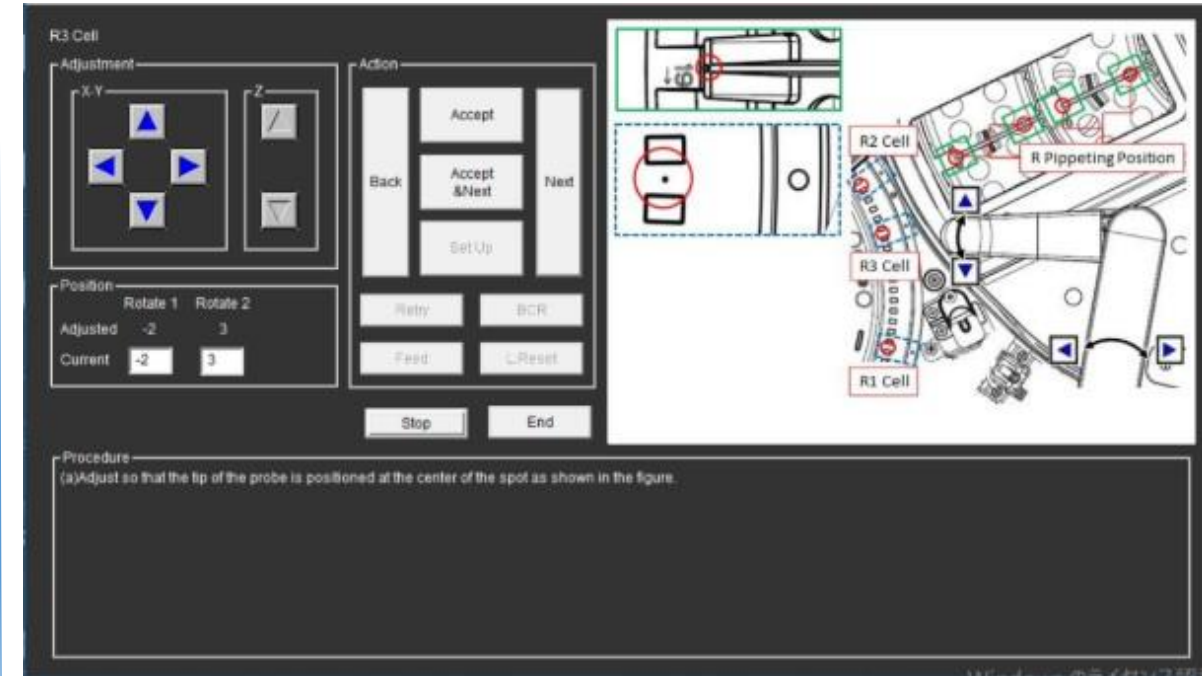


Position 6 (10/14)

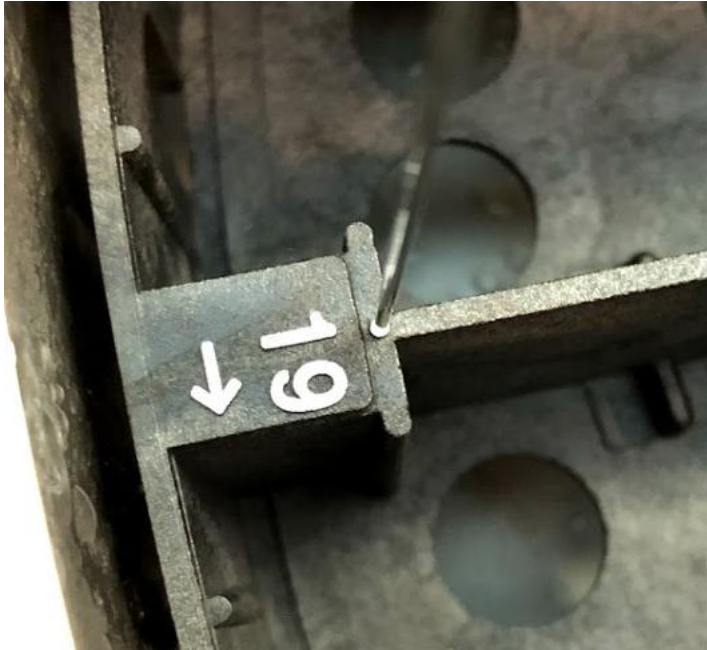


1. Adjust the probe , so that the tip of the probe is centered on the position.
2. Press 'Accept & Next' after adjustment.

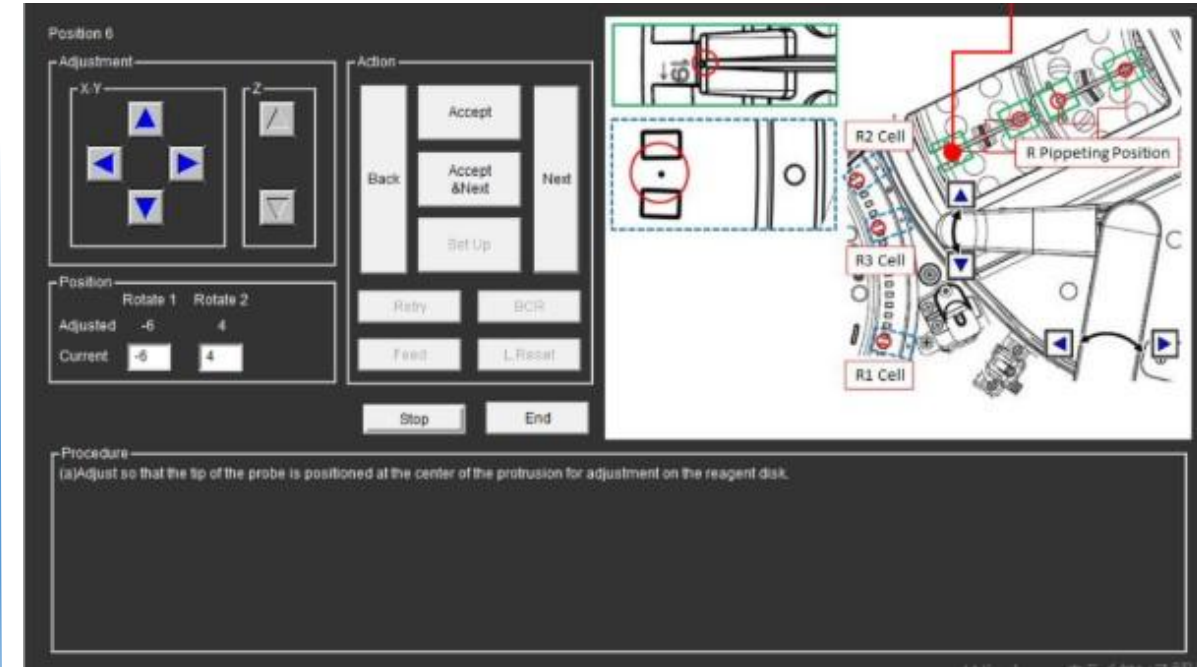

R Disk CW (11/14)

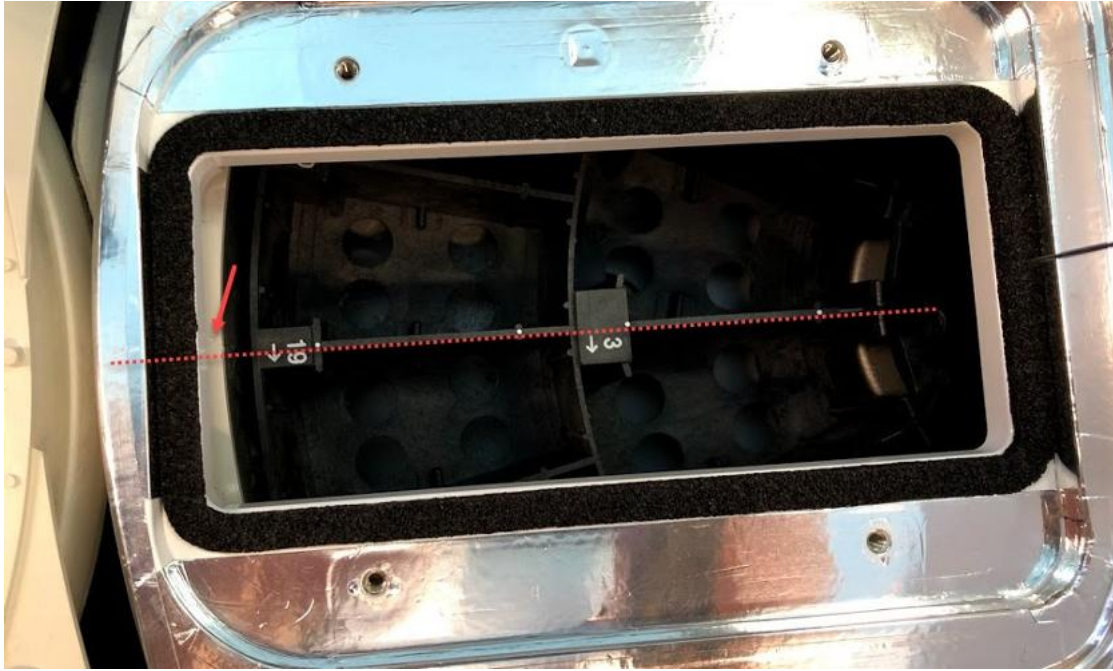


1. Adjust the Reagent Disk position , so that the Tip of the Probe is on the DOT.
2. Press 'Accept & Next' after adjustment.

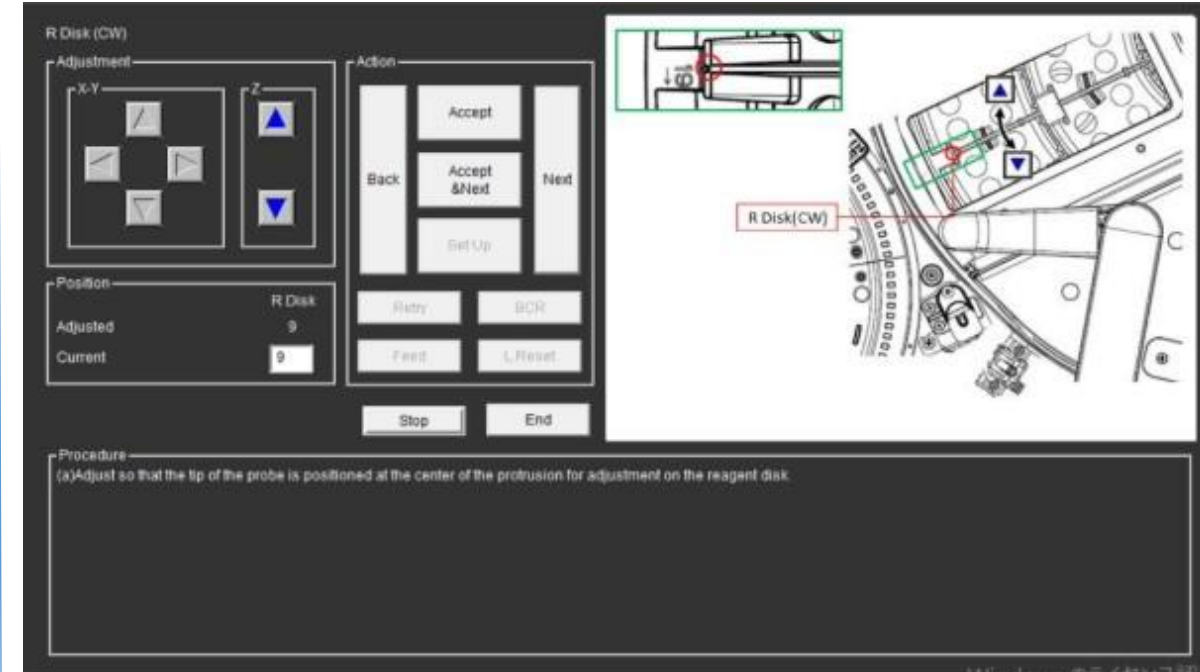


Pierce Wear Position(12/14)




1. Adjust the probe , so that it is centered as shown with the lines.
2. Press 'Accept & Next' after adjustment.

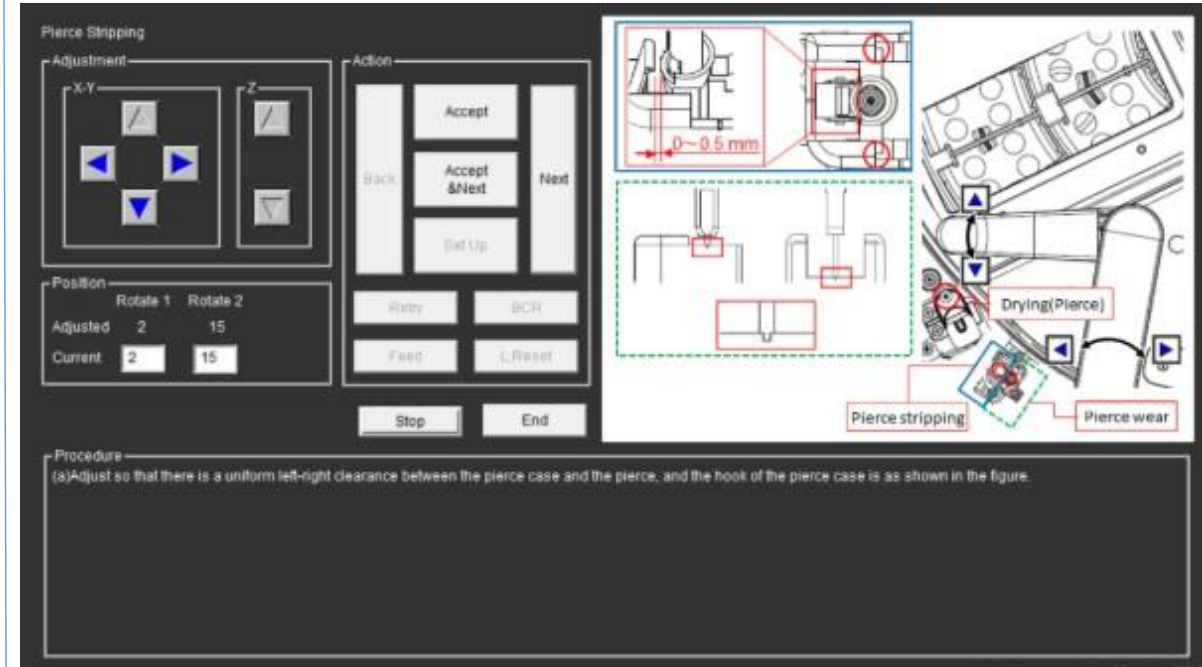

Drying (Pierce) (13/14)

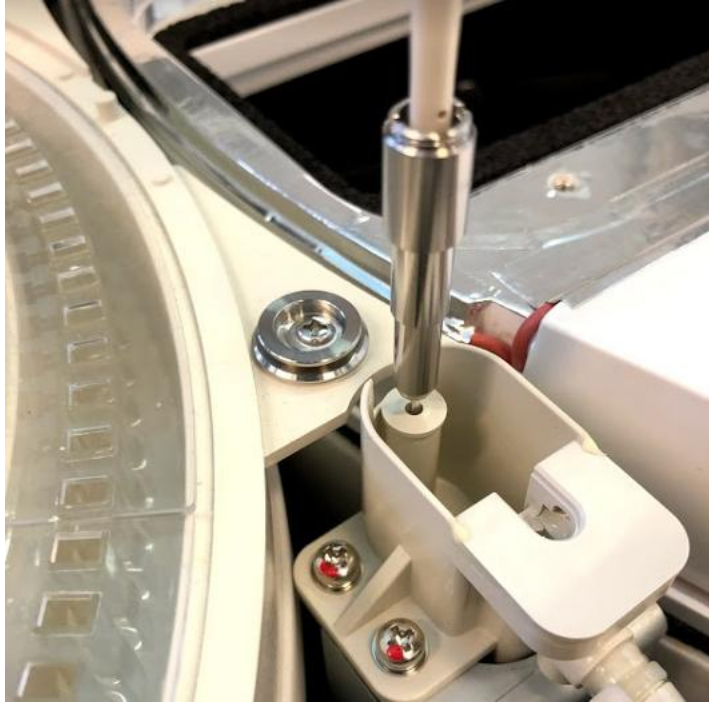


1. Adjust the piercer, so that it is centered at the drying hole.
2. Press 'Accept & Next' after adjustment.

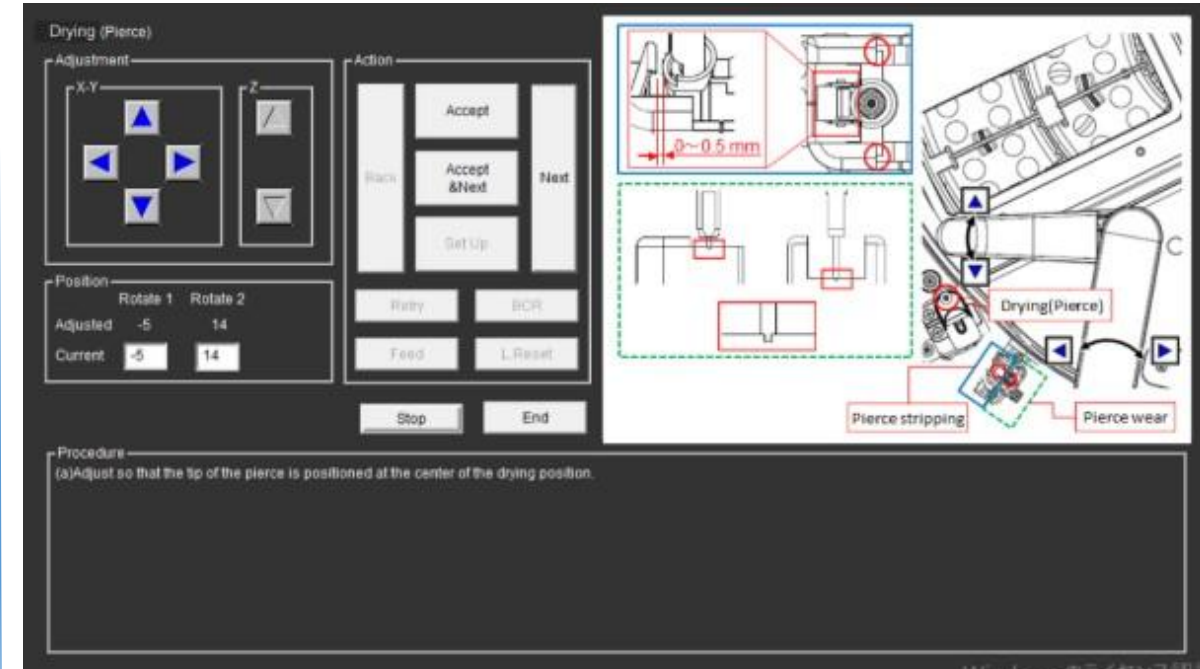

Pierce Stripping (14/14)



1. Adjust the piercer as shown above.
2. Press 'Accept' after adjustment.
3. Press 'End' to finish.

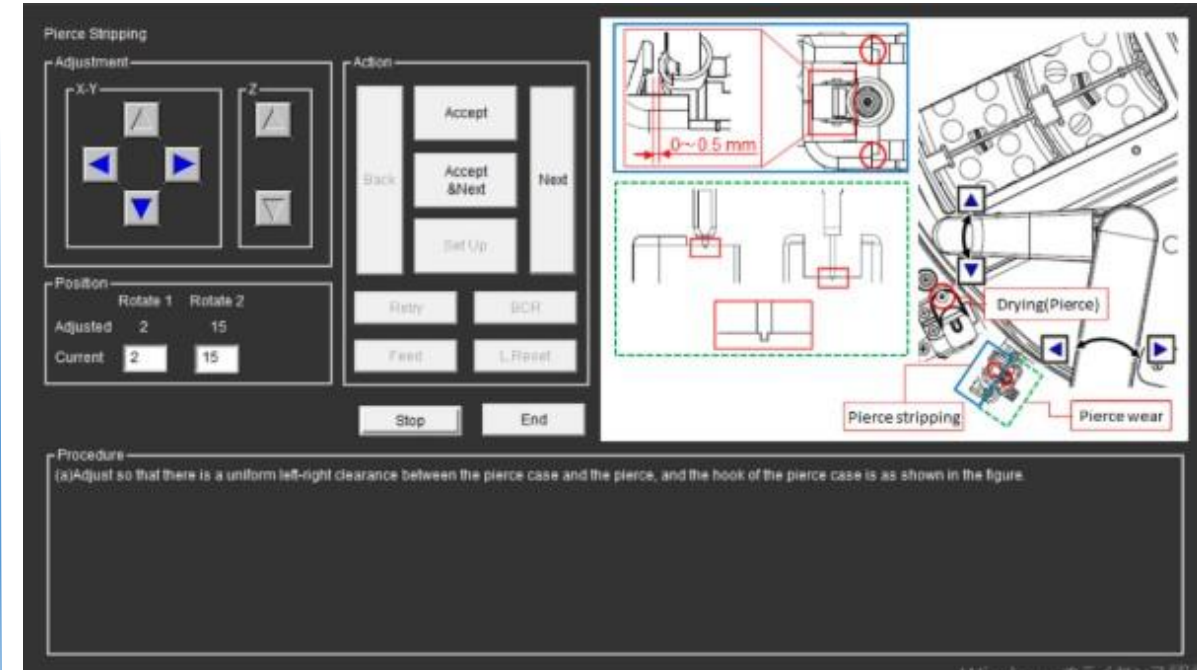

Reagent Probe Cell Vertical

Remove the cover shown in Figure 47 in advance.



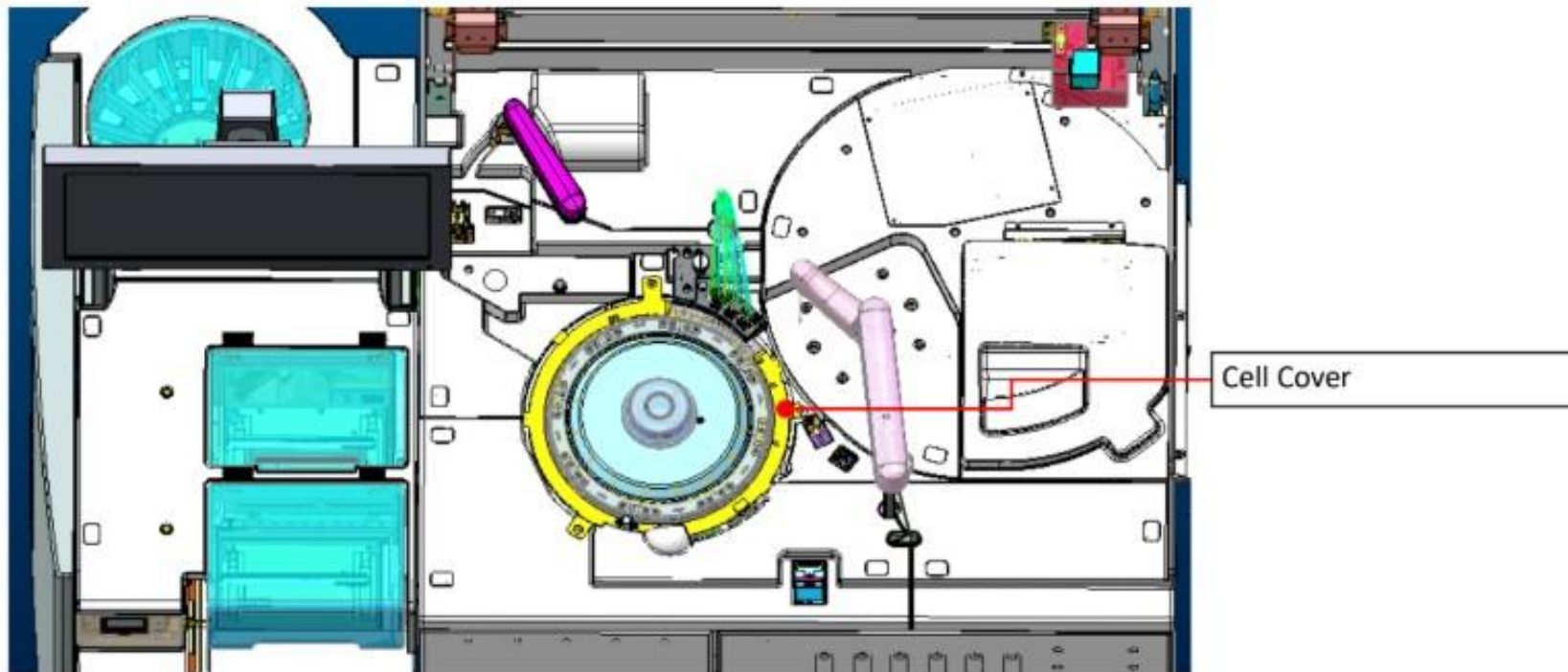


Figure 47. Removing the Cover

Cell Height (1/1)

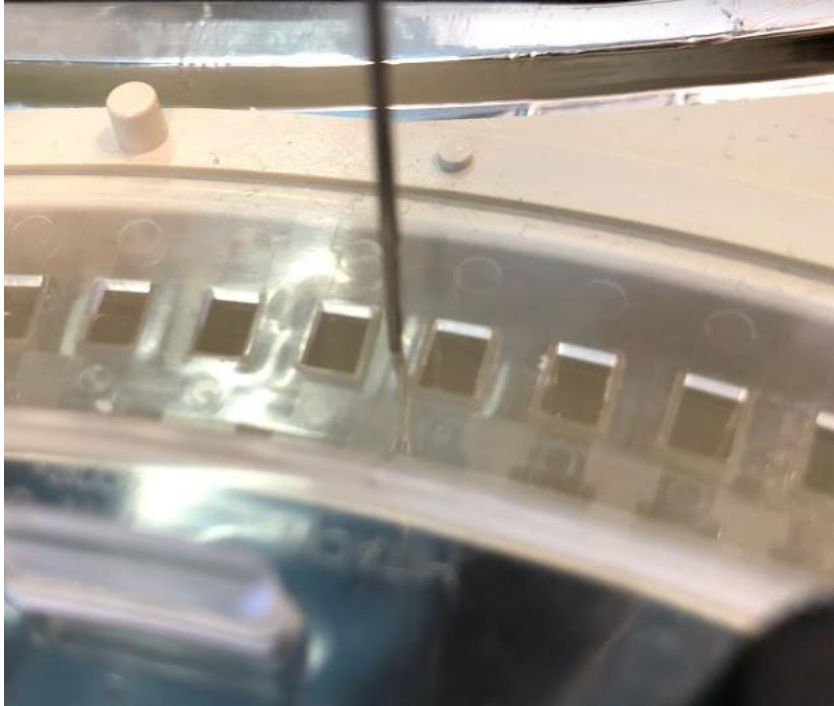

1. Press 'Setup' . Probe will move to adjustment position .
2. Adjust the tip of the probe , so that it slightly touches the cell surface
3. Press 'Accept' after adjustment.
4. Press 'End' to finish.

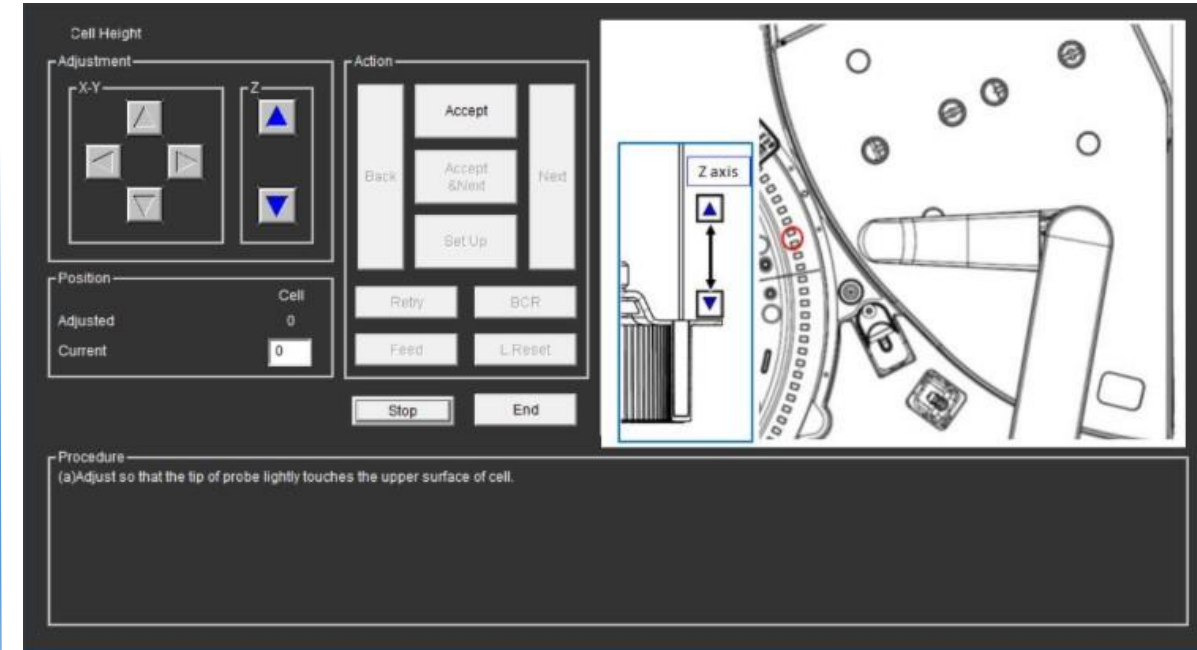

Reagent Probe Bottle Vertical


Prepare the adjustment Tool (3)

	Tool (3)
P/N	792-6449
Item name	Bottle Place C Set Assy
Appearance	 P/N 792-6453x2 Height adjustment position

Bottle Bottom (1/1) – Preparation :



1. Open the Reagent Disk cover and place adjustment jigs to position 14 & 40.
2. Press 'Setup'

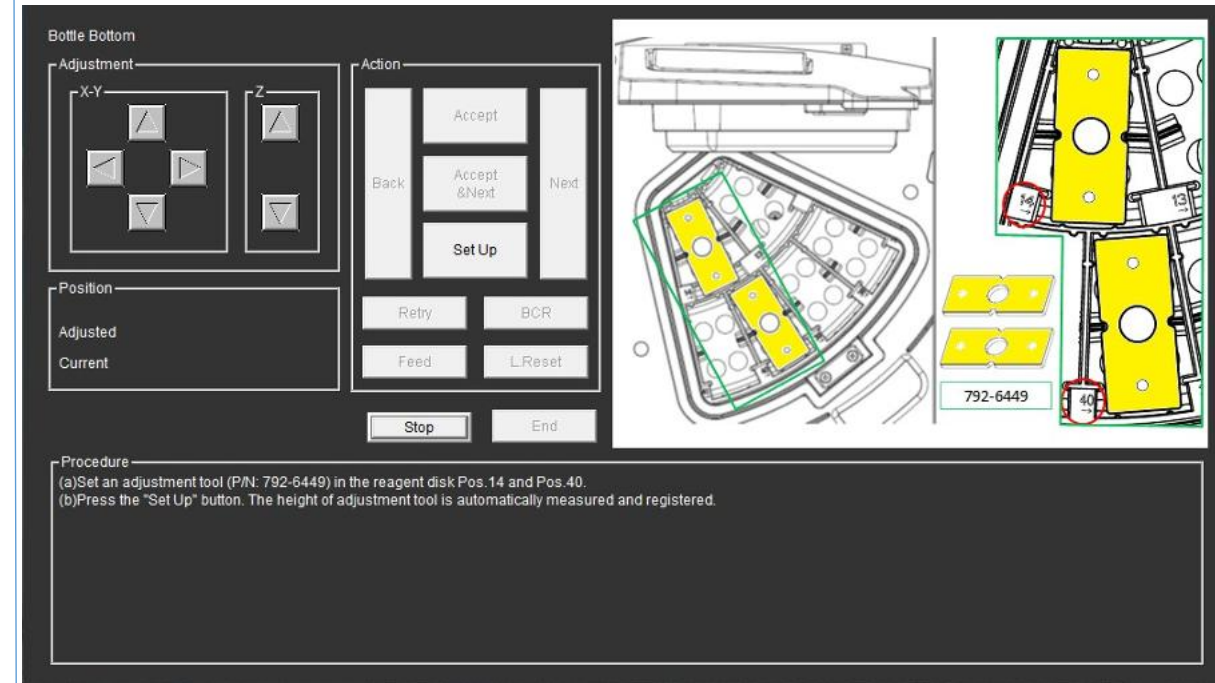


Adjustment Flowchart – Start Page

Bottle Bottom (1/1)




1. Probe will move automatically , first to position 14 , then to position 40 .

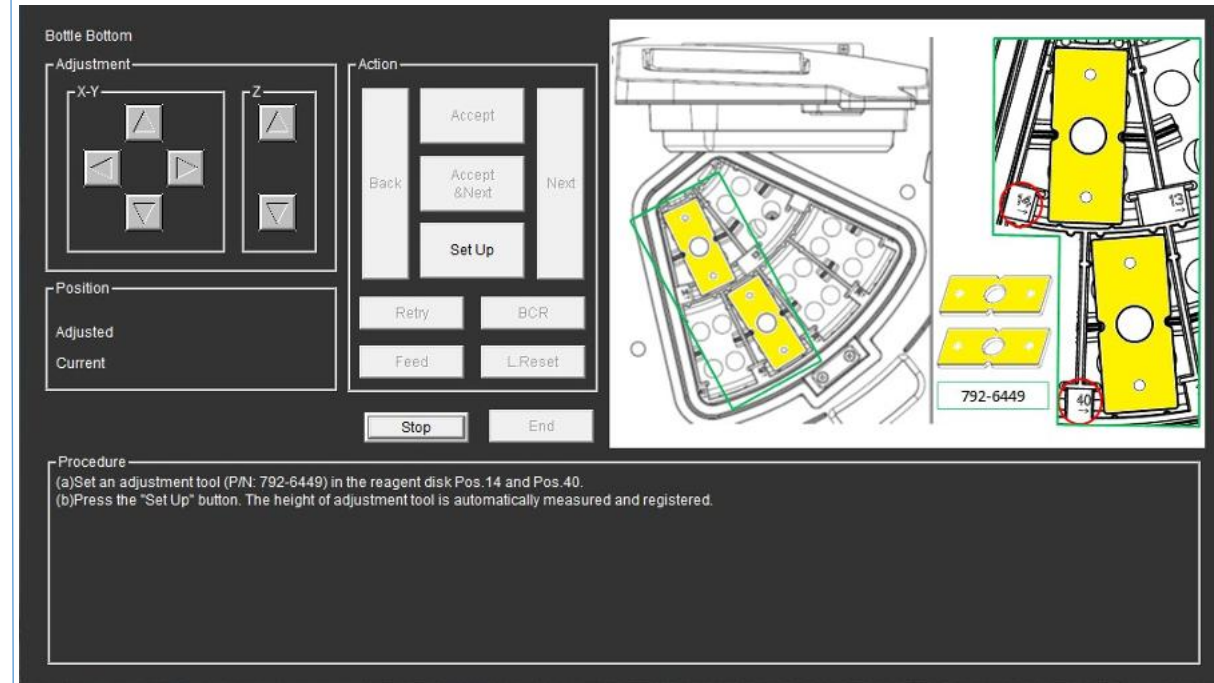


Bottle Bottom (1/1) :



1. Remove the adjustment Jigs.
2. Press 'Setup'.
3. Press 'End' to finish

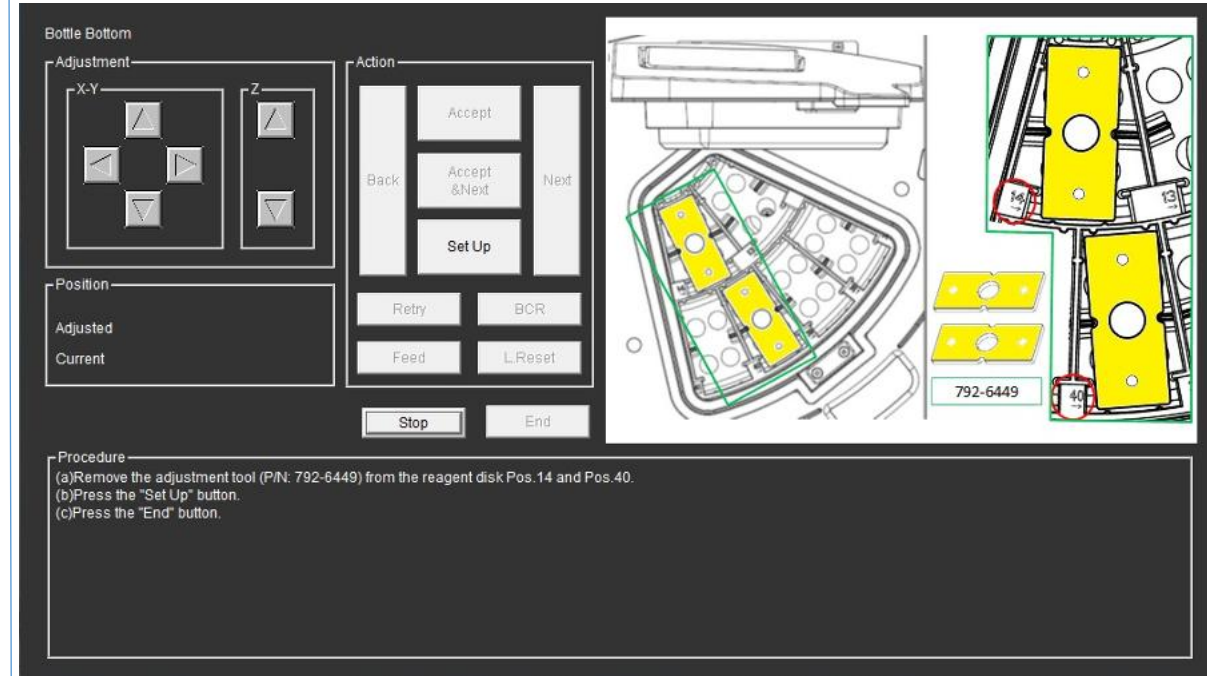

Reaction Rinse Vertical

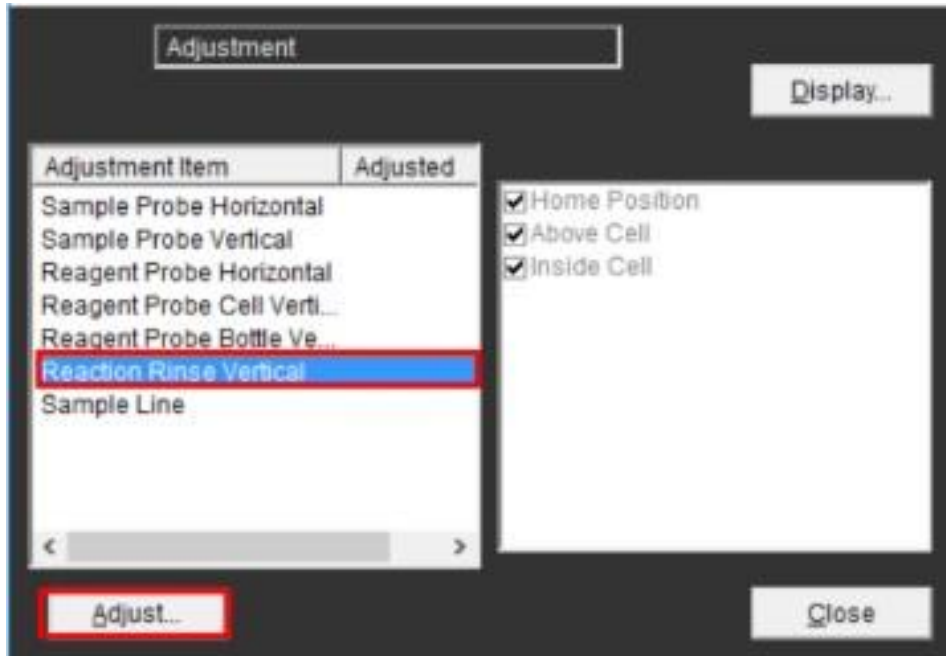


Rinse Mechanism will stop above 1mm from the reaction cell surface. Then every press to 'Next' button is the next adjustment position.

Home Position (1/3)




Check the position of the probes.

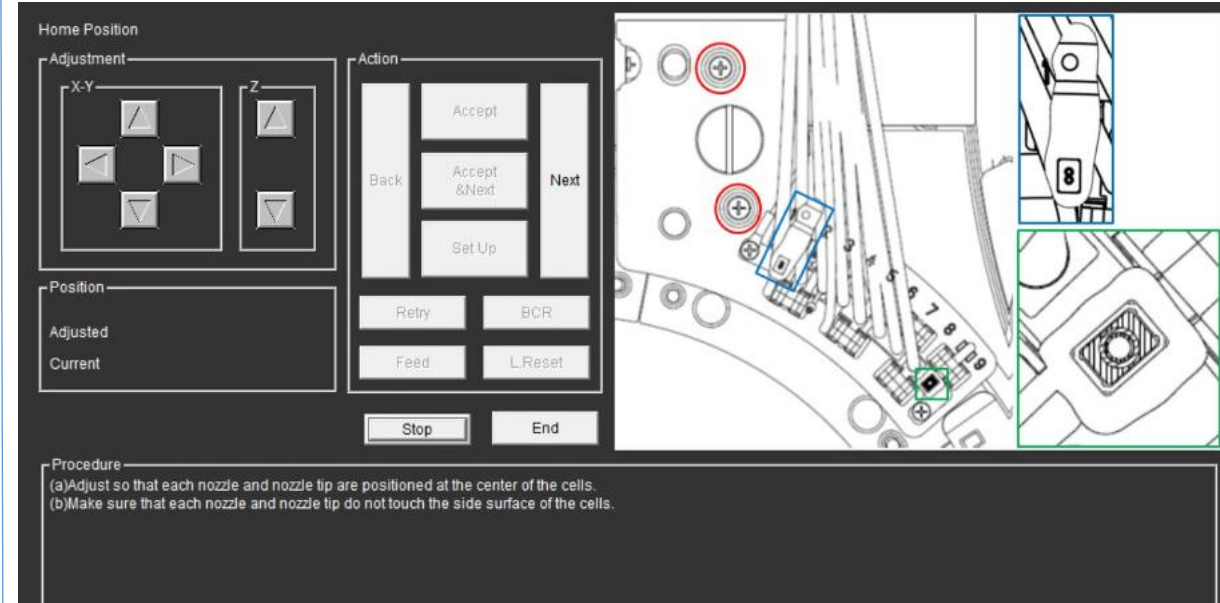

Cell Top (2/3)





Check the position of the probes.

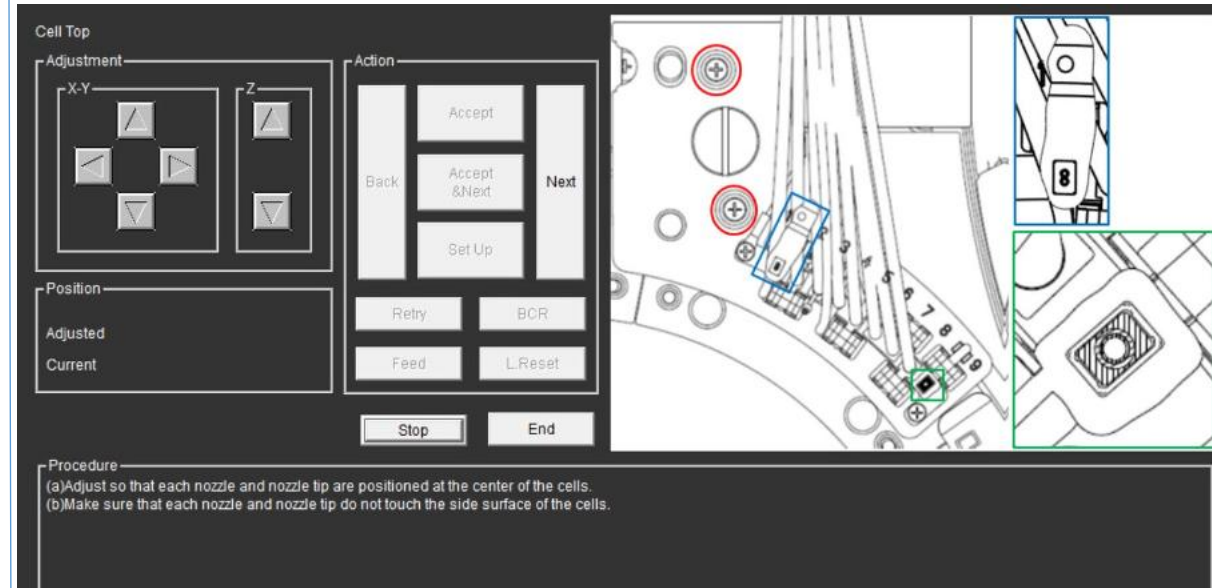

Over Flow (3/3)




Check the position of the probes. If needed , adjust by the fixing screws shown.

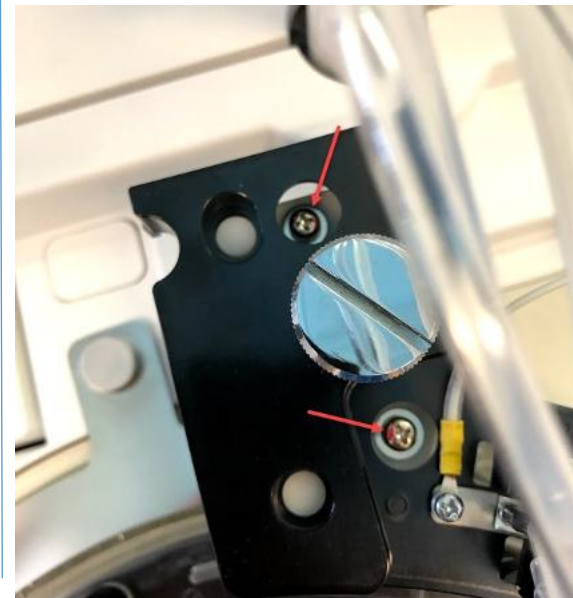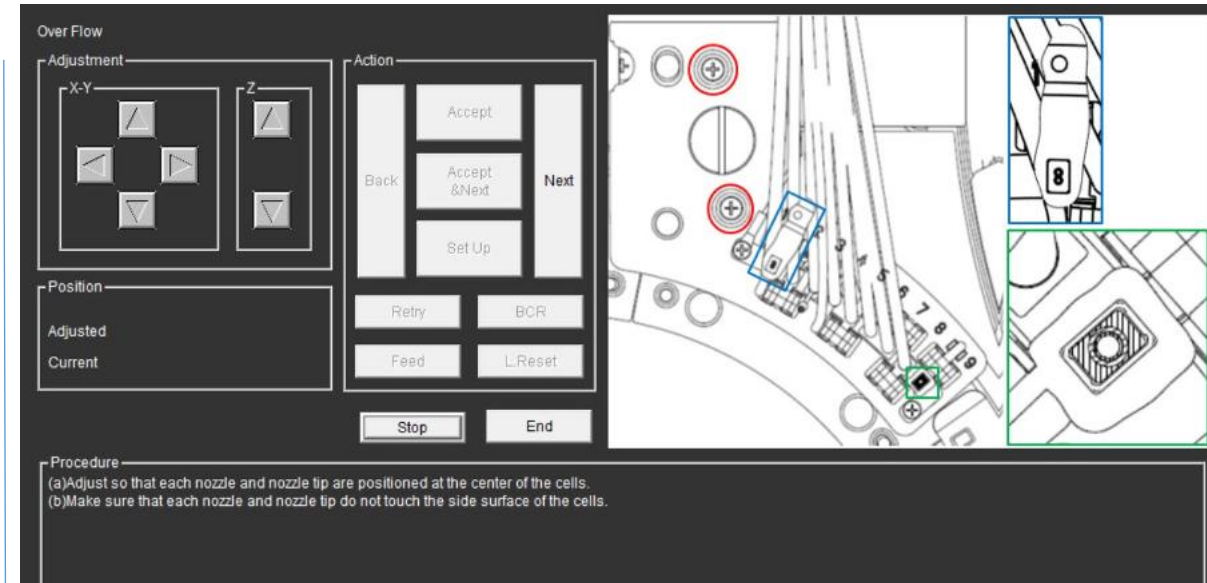

[Adjustment Flowchart – Start Page](#)

Sample Line

Prepare the adjustment Tool (4)

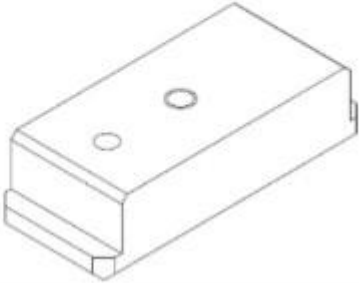	Tool(4)
P/N	837-6416
Item name	Rack Rotor ADJ Jig
Appearance	

Remove the cover shown in Figure 54 in advance.


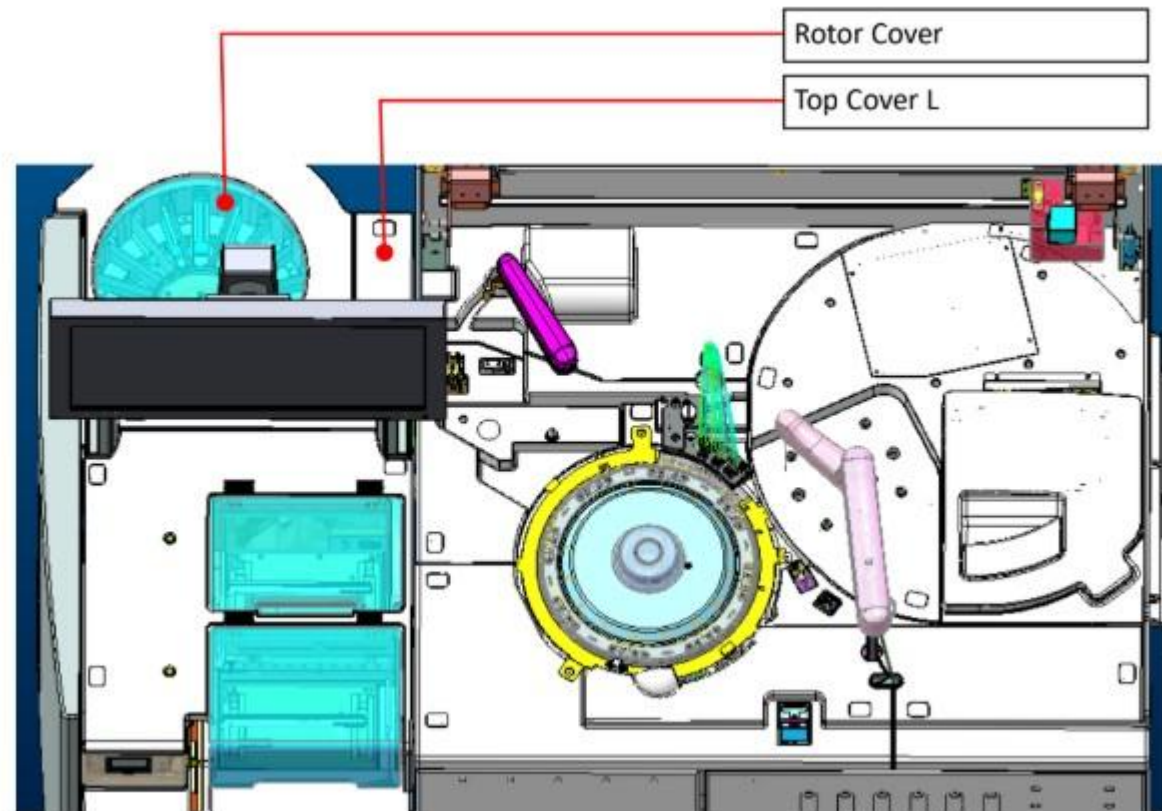


Figure 54. Removing the Cover

Rack Transfer End (1/2)

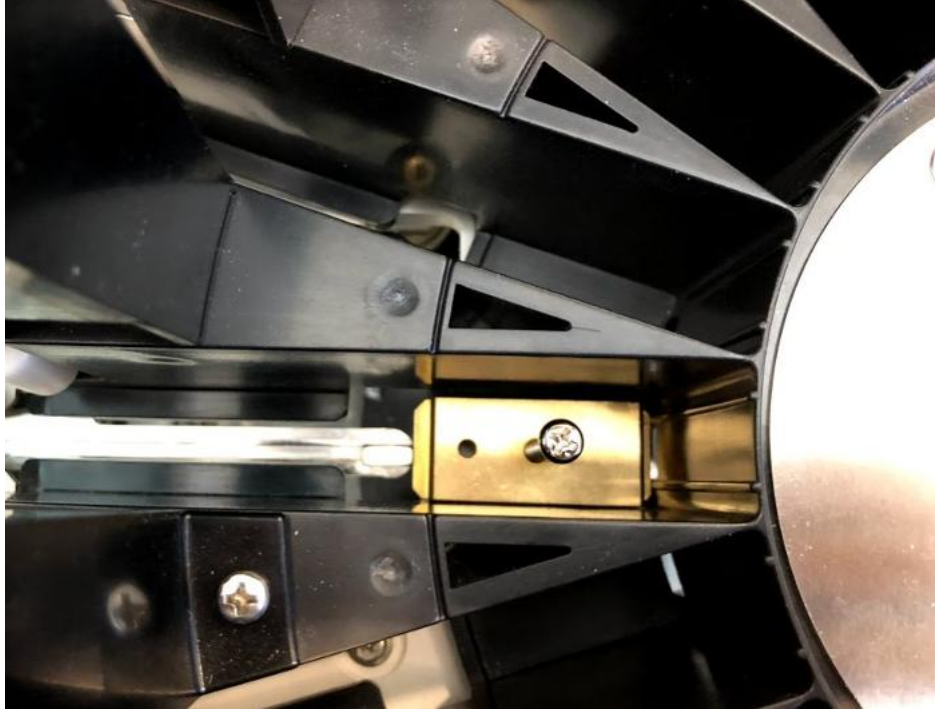

1. Press 'Set Up' to move Pin to rotor.
2. Set the Jig, so that it is slightly touching Pin.
3. Remove the Jig.
4. Press 'Accept & Next' after adjustment.

Rack Transfer End

Adjustment

X-Y

Z

Position

Adjusted -6

Current -6

Action

Back

Accept

Accept & Next

Next

Set Up

Retry

BCR

Feed

L Reset

Stop

End


Procedure

(a) Press the "Set Up" button to move the hook of the sample line toward the rack rotor.

(b) Set an adjustment tool (P/N: 837-6415) as shown in the figure.

(c) Adjust so that there is no gap between the hook of the sample line and the adjustment tool.

(d) After adjustment, remove the adjustment tool from the rack rotor.

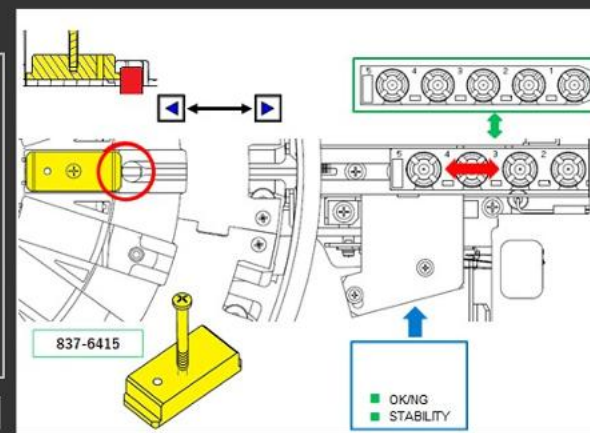

(BCR Adjustment) Waiting Position (2/2)



Figure 58. Adjustment of BCR

LED

Mounting position
(2 screws)



1. Place a Sample Rack on setting position.
2. Press 'BCR' and check the status of two LED's , both needs to be Green.
3. Remove the Rack.
4. Press 'Accept' after adjustment.
5. Press 'End' to finish.

BCR

Adjustment

X-Y

Z

Position

Adjusted

Current

Action

Accept

Accept & Next

Set Up

Back

Next

Retry

BCR

Feed

LReset

Stop

End

Procedure

(a) Press the "Next" button to move the hook of the sample line to waiting position.

(b) Set the rack for f013 test tube on the waiting position of the sample line.

(c) Press the "BCR" button to turn on the BCR.

(d) Move the rack slightly left and right by the clearance between the hook of the sample line. While moving the rack, make sure that all LEDs are lit in green. If it is not so, adjust the mounting position of the BCR.

(e) Remove the rack for f013 test tube. Set the rack for f016 test tube.

(f) Make sure that all LEDs are lit in green, just as in the case of the rack for f013 test tube. If the mounting position of the BCR is changed for adjustment of the rack for f016 test tube, check the rack for f013 test tube again.

(g) Press the "BCR" button to turn off the BCR.

(h) Press the "End" button.

Adjustment Flowchart – Start Page